
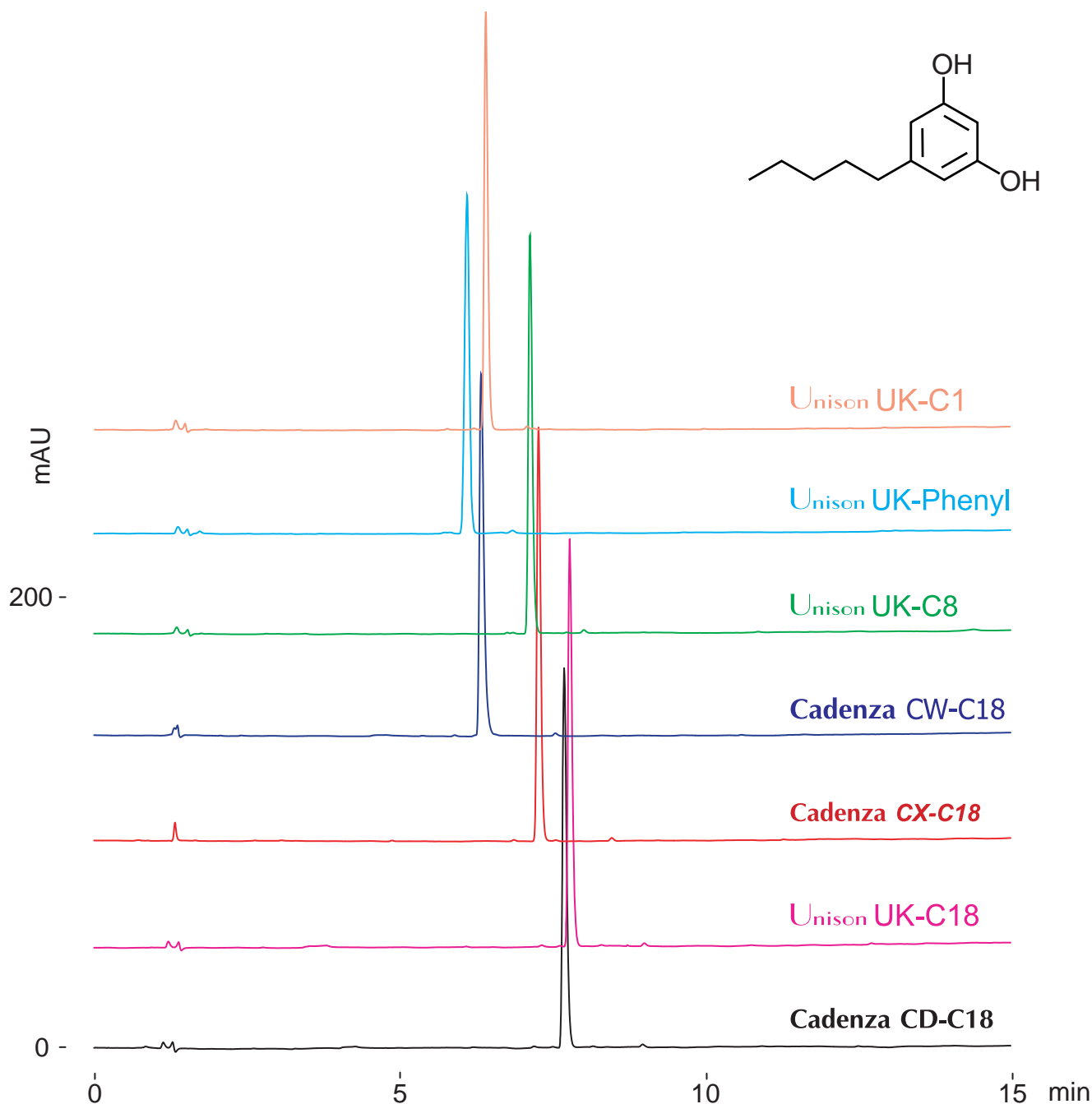


Unison UK-C1
 Unison UK-Phenyl
 Unison UK-C8
 Unison UK-C18
 Cadenza CX-C18
 Cadenza CW-C18
 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Olivetol (Phenols)
 オリベトール (フェノール類)

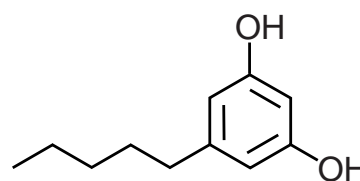

100 x 4.6 mm
 A: water / acetonitrile / formic acid = 90 / 10 / 0.1
 B: water / acetonitrile / formic acid = 10 / 90 / 0.1
 15-100%B (0-15min)
 1.0mL/min (11-16MPa), 37°C, 260nm
 10µL (1.0mg/mL, water / acetonitrile / formic acid = 90 / 10 / 1)