
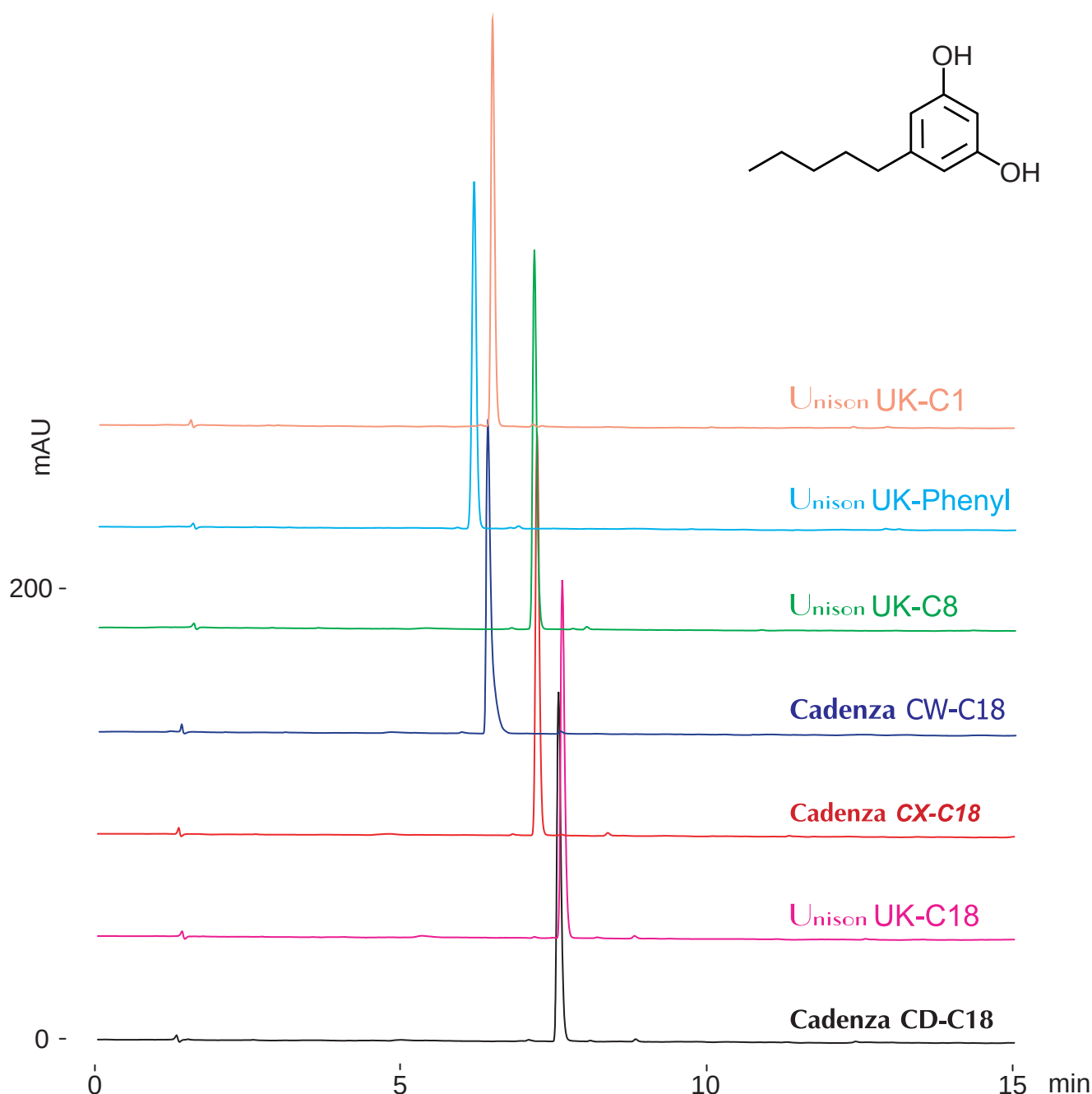


Unison UK-C1  
Unison UK-Phenyl  
Unison UK-C8  
Unison UK-C18  
Cadenza CX-C18  
Cadenza CW-C18  
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Olivetol (Phenols) オリベトール (フェノール類)

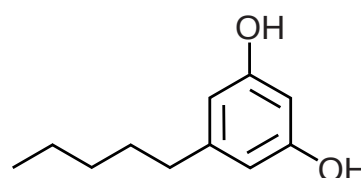

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

15-100%B (0-15min)

1.0mL/min (11-15MPa), 37°C, 260nm

10μL (1.0mg/mL, 100mM ammonium formate / acetonitrile = 90 / 10)