

Unison UK-Amino


Nardis ND-NX

150 x 4.6 mm

HPLC COLUMN APPLICATION

2'-deoxy-2'-fluoroguanosine, 2'-deoxyguanosine hydrate

2'-デオキシ-2'-フルオログアノシン、2'-デオキシグアノシン水和物

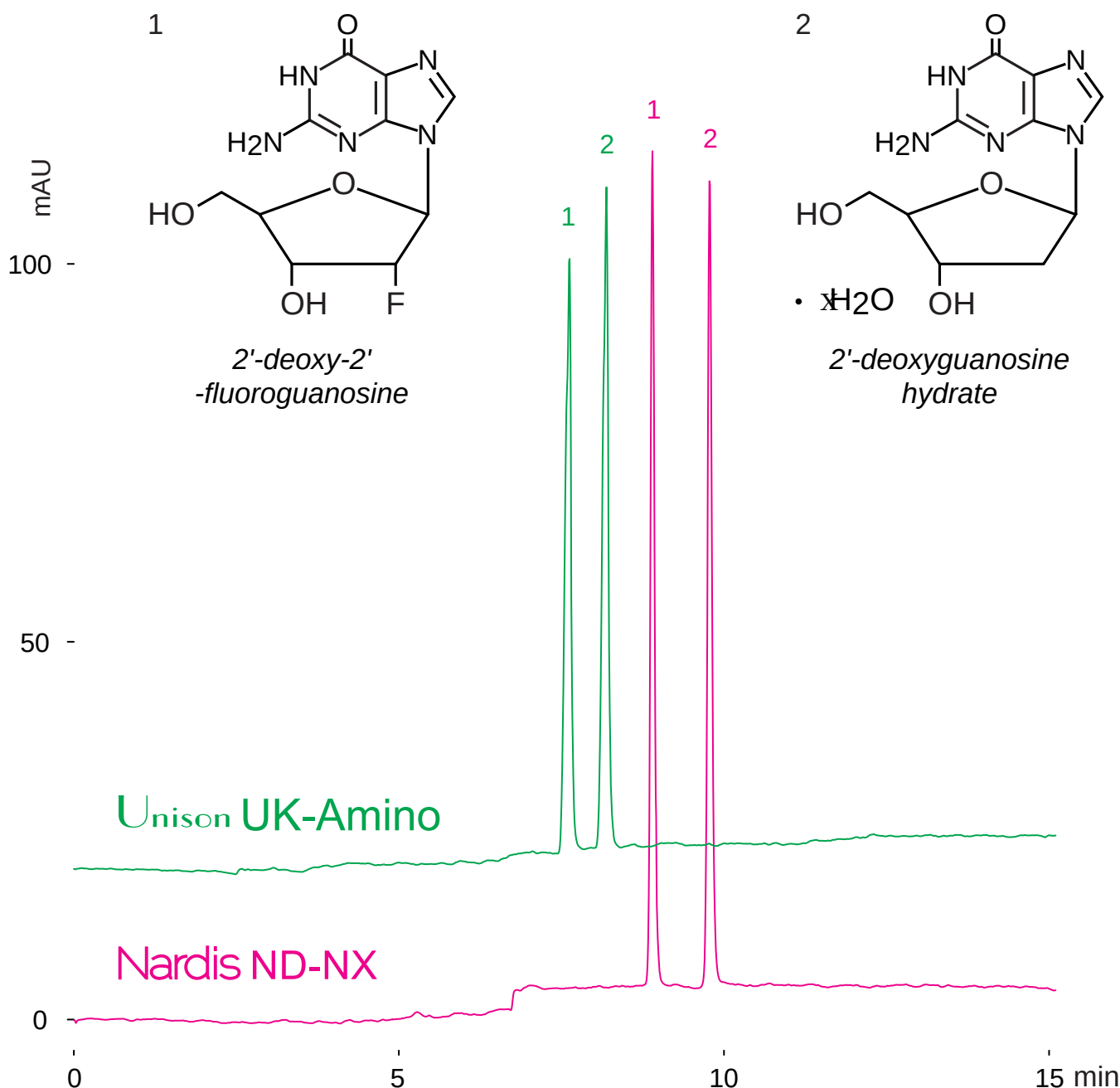

150 x 4.6 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium acetate = 50 / 50

0-100%B (0-15min)

1.0mL/min (6MPa), 37°C, 260nm

10μL (50μg/mL, 0.1M ammonia solution)