


Nardis ND-RX
Scherzo SW-C18
Scherzo SM-C18
Scherzo SS-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Chloroquine analogues (antimalarial)


クロロキン類縁体 (抗マラリア薬)

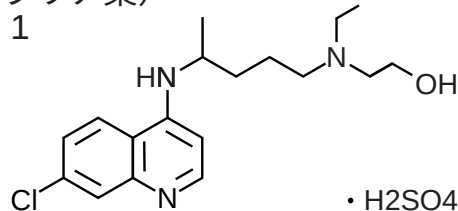

hydroxychloroquine sulfate(0.1mg/mL)

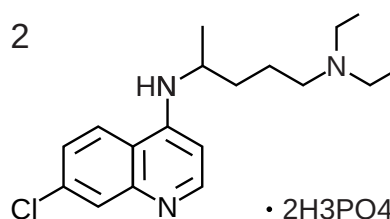

chloroquine diphosphate(0.1mg/mL)

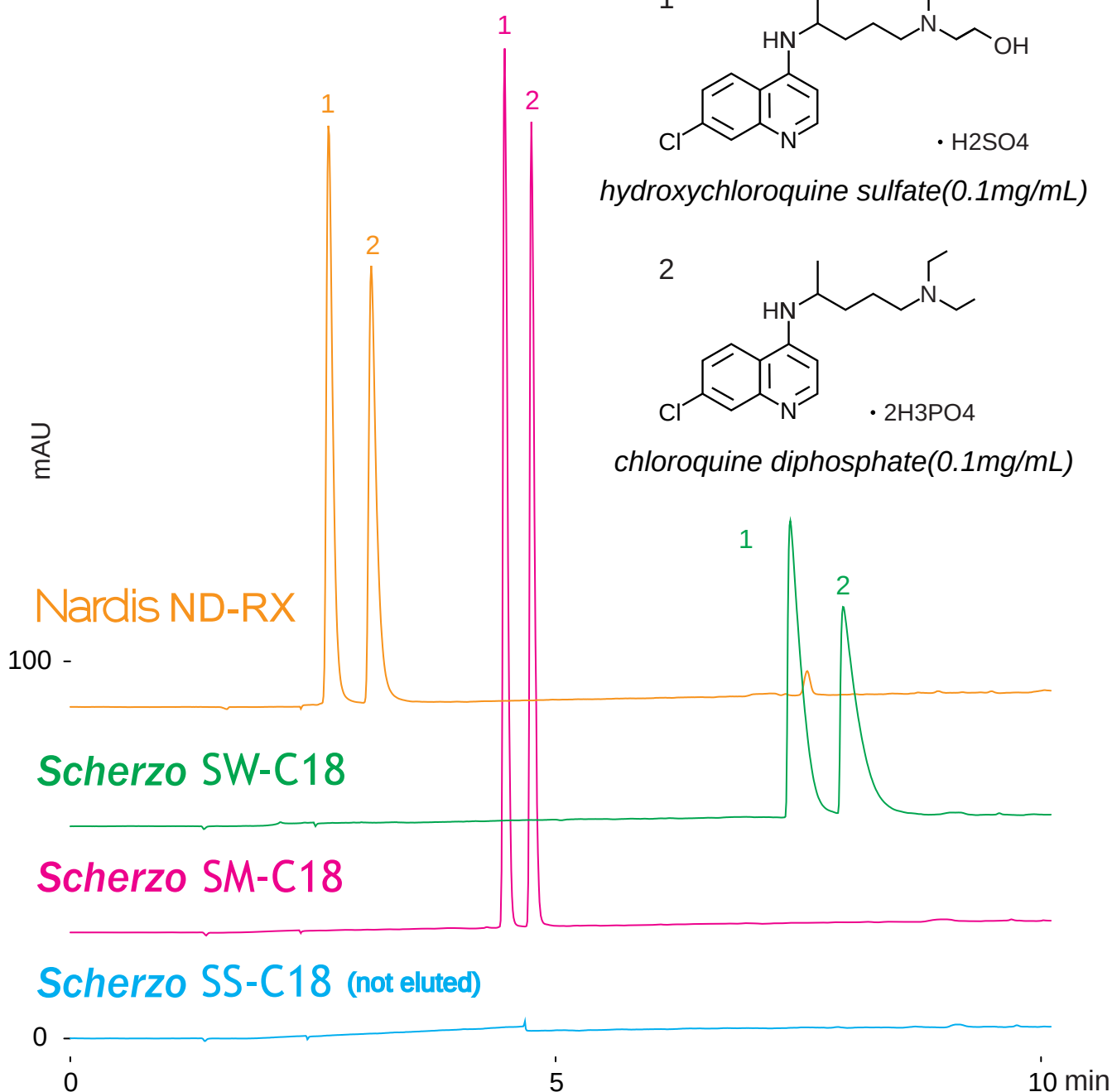

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-100%B (0-10min)

1.0mL/min (7-8MPa), 37°C, 260nm

10μL (0.1mg/mL, water)