

Nardis ND-NX MF


Nardis ND-NX

150 x 2 mm

HPLC COLUMN APPLICATION

Hydroxypicolinic acid isomers

ヒドロキシピコリン酸異性体

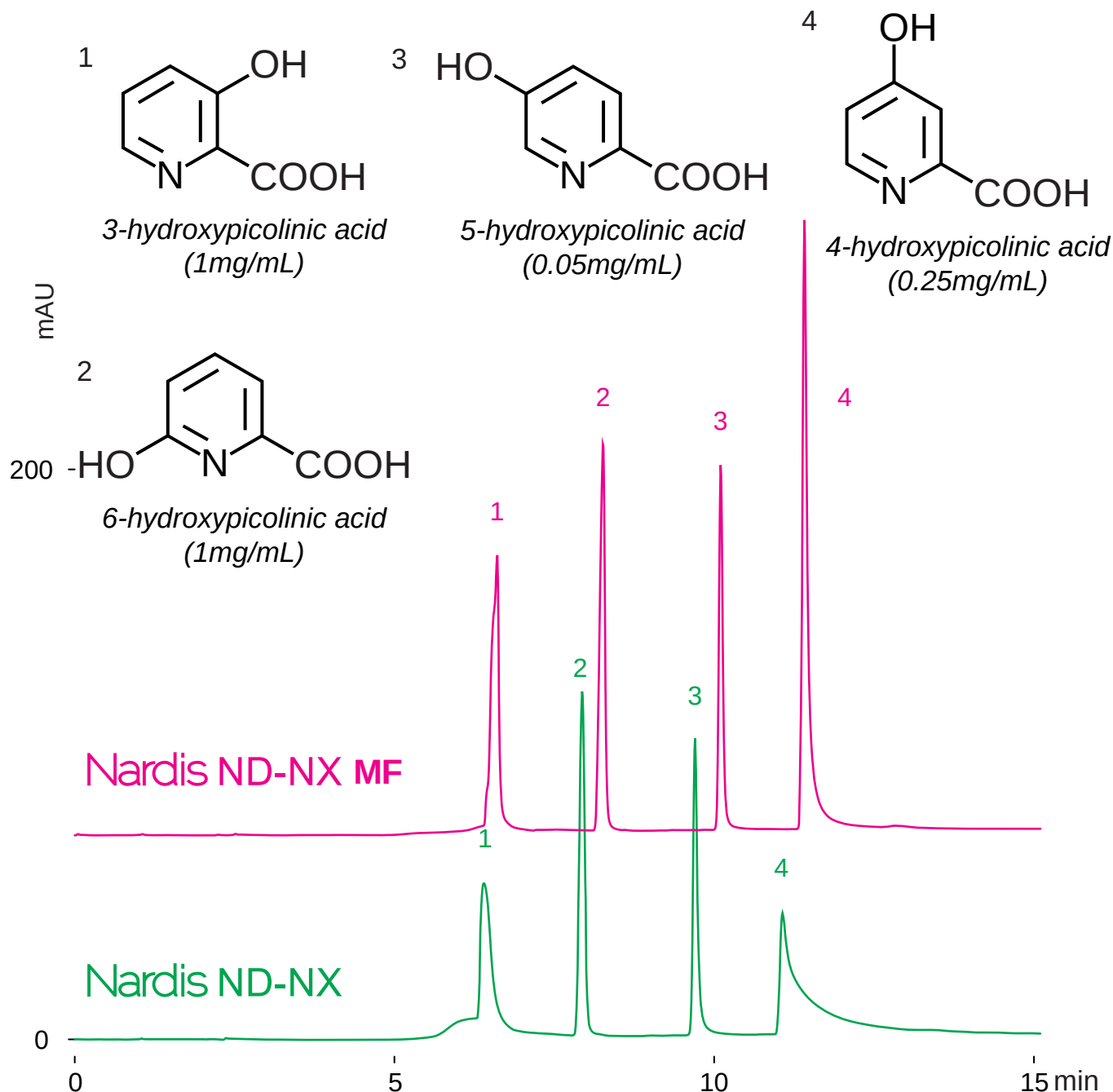

150 x 2 mm

A: acetonitrile / 100mM ammonium acetate = 95 / 5

B: acetonitrile / 100mM ammonium acetate = 50 / 50

0-70%B (0-15min)

0.2mL/min (5-9MPa), 37°C, 260nm

2μL (1-0.05mg/mL, ACN/100mM AcONH4=50/50)