

Cadenza CW-C18

Cadenza CX-C18 Unison UK-C1

Unison UK-C18 Unison UK-Phenyl


Cadenza CD-C18 Unison UK-C8

100 x 4.6 mm


HPLC COLUMN APPLICATION

thiazoyl blue tetrazolium bromide (MTT assay)

臭化チアゾイルブルーテトラゾリウム (MTTアッセイ)

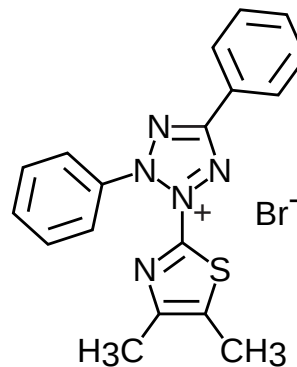

thiazoyl blue tetrazolium bromide

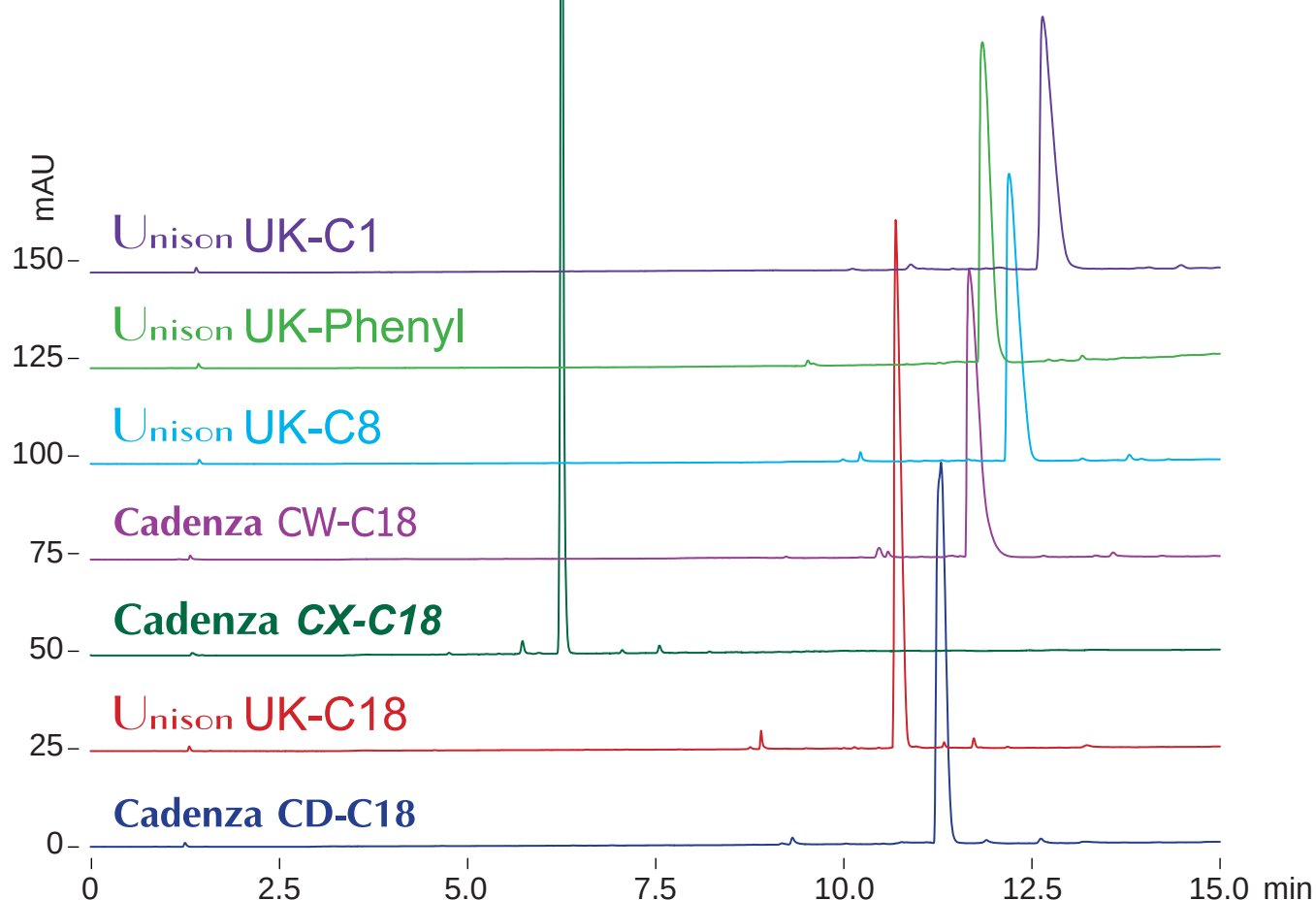

100 x 4.6 mm

A: acetonitrile / water / formic acid = 5 / 95 / 0.1

B: acetonitrile / water / formic acid = 95 / 5 / 0.1

0%B (0-1min), 0-60%B (1-15min)

1mL/min (10-12MPa), 37°C, 260nm

5μL (100μg/mL, water / formic acid = 100 / 1)