


Cadenza CX-C18

150 x 3 mm

HPLC COLUMN APPLICATION

thiazoyl blue tetrazolium bromide (MTT assay)

臭化チアゾイルブルーテトラゾリウム (MTTアッセイ)

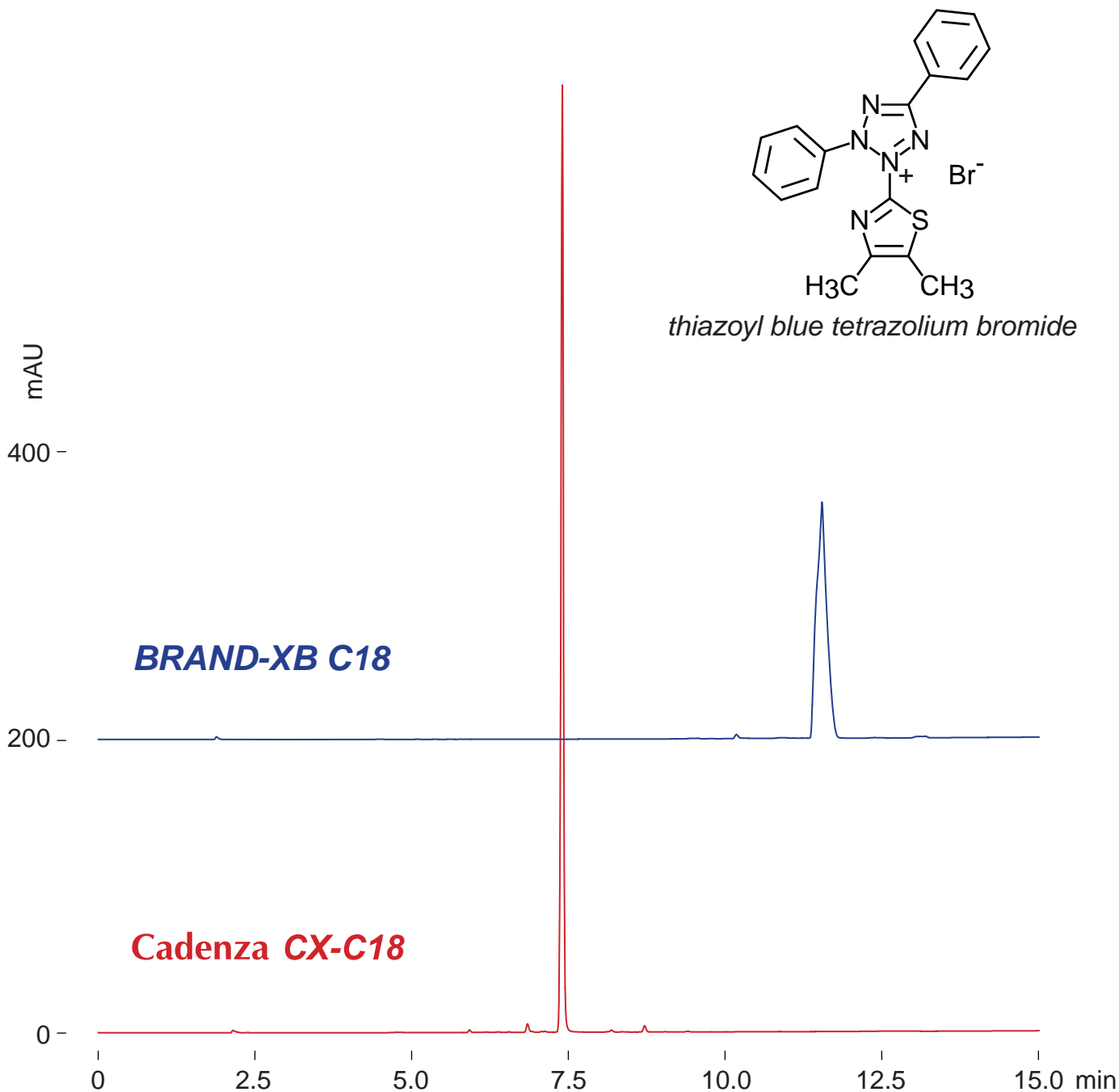

150 x 3 mm

A: acetonitrile / water / formic acid = 5 / 95 / 0.1

B: acetonitrile / water / formic acid = 95 / 5 / 0.1

0%B (0-1min), 0-60%B (1-15min)

0.4mL/min (10MPa), 37°C, 260nm

5µL (100µg/mL, water / formic acid = 100 / 1)