

Cadenza CW-C18

Cadenza CX-C18 Unison UK-C1

Unison UK-C18 Unison UK-Phenyl


Cadenza CD-C18 Unison UK-C8

100 x 4.6 mm

HPLC COLUMN APPLICATION

thiazoyl blue tetrazolium bromide (MTT assay)

臭化チアゾイルブルーテトラゾリウム (MTTアッセイ)

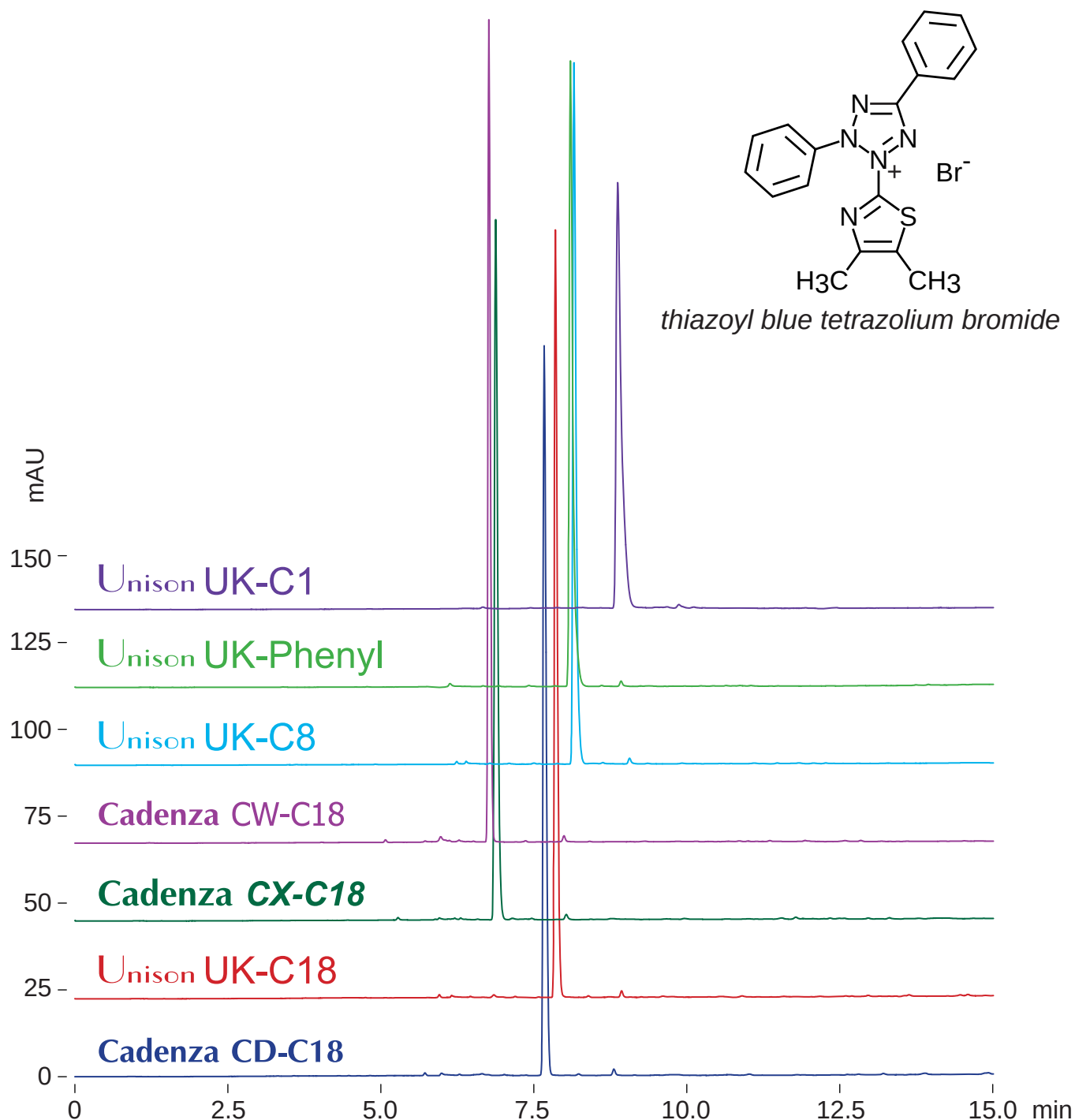

100 x 4.6 mm

A: acetonitrile / 100mM ammonium formate = 5 / 95

B: acetonitrile / 100mM ammonium formate = 95 / 5

10-100%B (0-15min)

1mL/min (10-12MPa), 37°C, 260nm

5μL (100μg/mL, 100mM ammonium formate)