

Scherzo SM-C18 Nardis ND-RX


Scherzo SW-C18 Scherzo SS-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

(S) - (+) -2- (Anilinomethyl) pyrrolidine(Proline-type organocatalyst)

(S) -(+) -2 -(アニリノメチル) ピロリジン (プロリン型有機触媒)

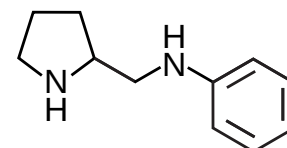

(S) - (+) -2- (anilinomethyl)
pyrrolidine

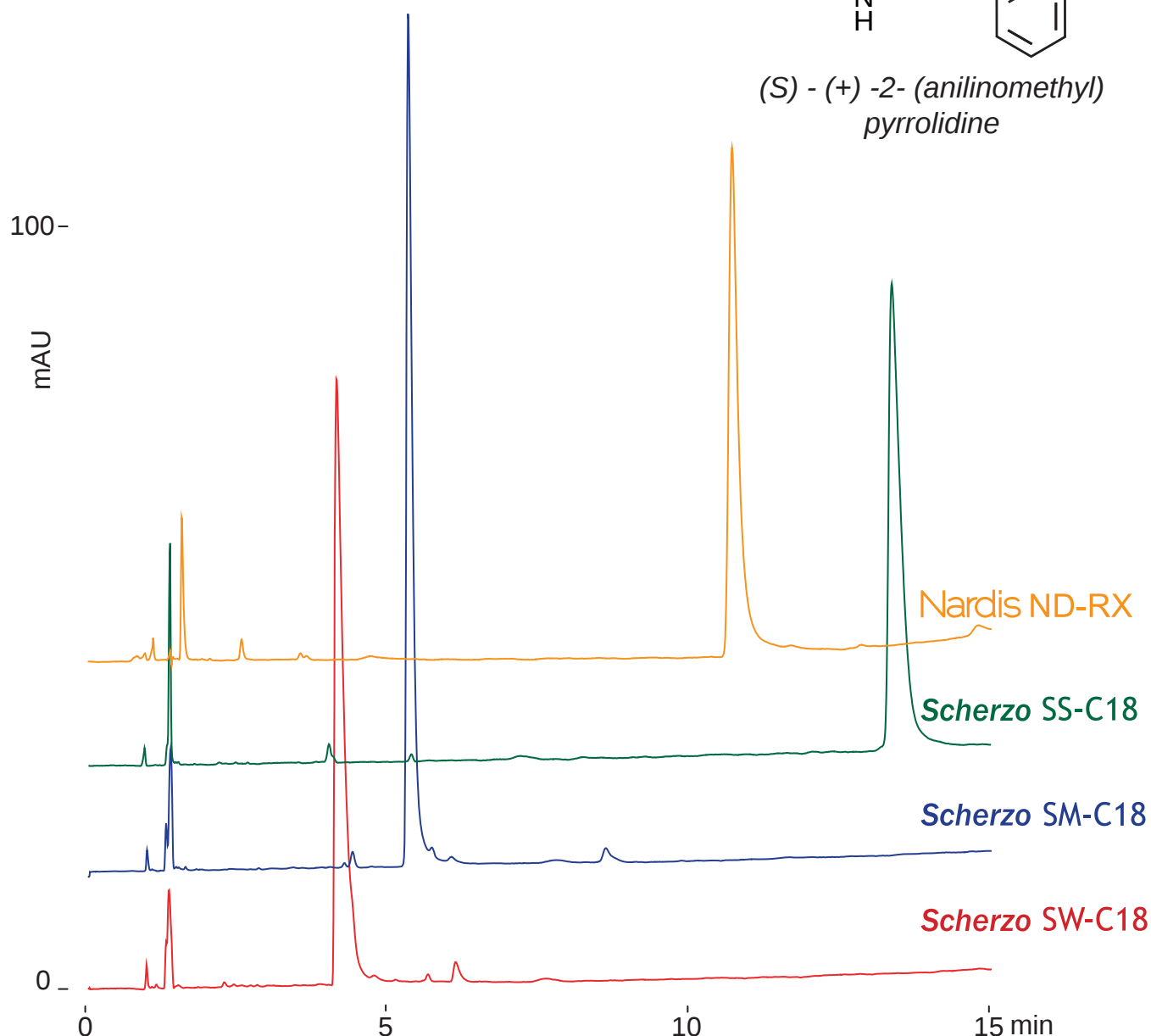

100 x 4.6 mm

A: 10mM ammonium formate / acetonitrile = 95 / 5

B: 100mM ammonium formate / acetonitrile = 50 / 50

30-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 245nm

10μL (0.05mg/mL, 100mM ammonium formate / acetonitrile = 90 / 10)