


Nardis ND-NX

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Guanidine hydrochloride salt

グアニジン塩酸塩の LC-MS 分析

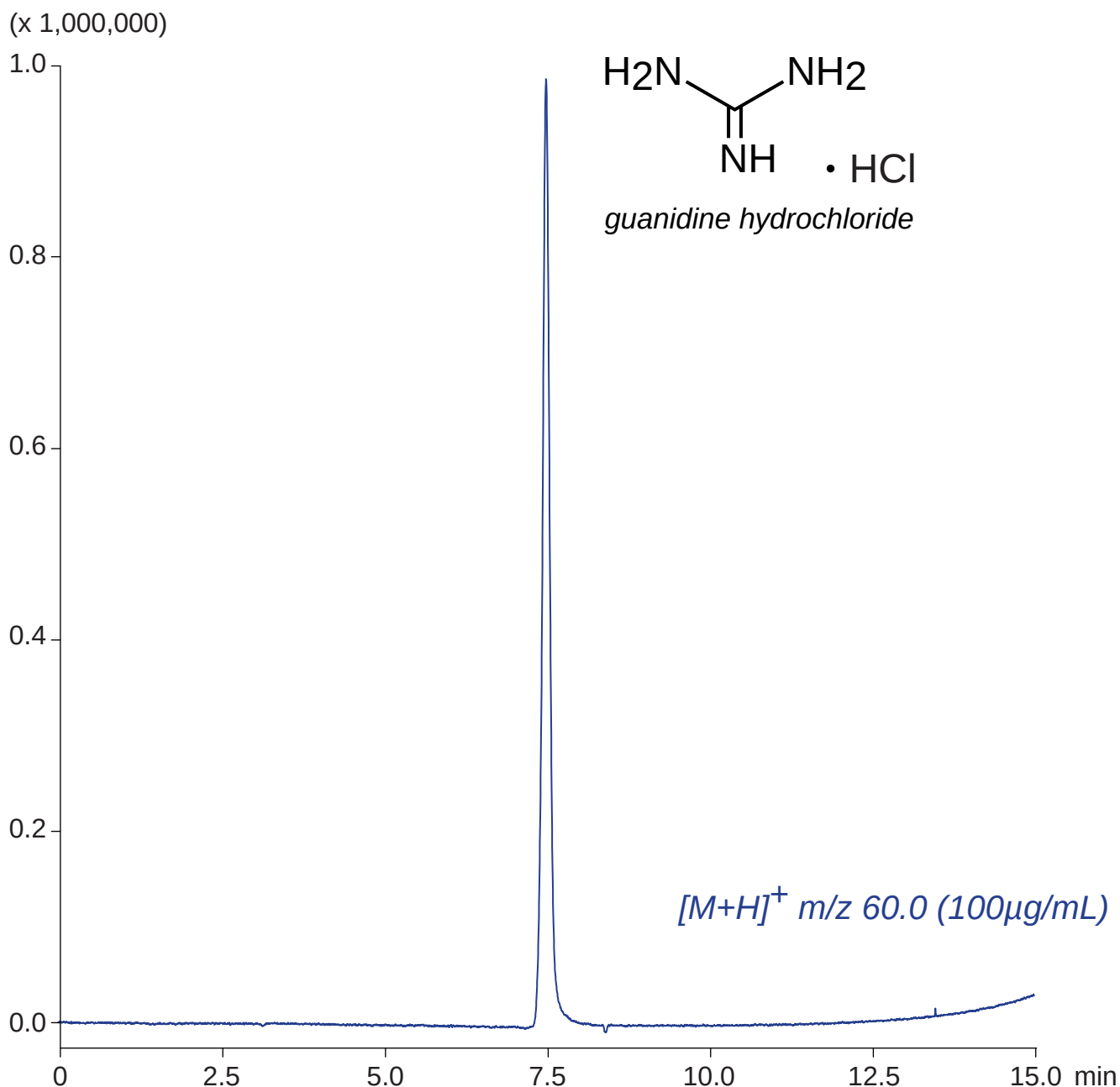

Nardis ND-NX, 150 x 3 mm

A: acetonitrile / 100mM ammonium formate = 90 / 10

B: acetonitrile / water / formic acid = 20 / 80 / 0.5

15-100%B (0-15min)

0.4mL/min (5MPa), 37°C, Single Quad. MS (ESI)

5 μL (50 μg/mL, 100mM AcONH₄ / ACN = 30 / 70)