


Unison UK-C1 Cadenza GL-C18
 Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

2- (4-Hydroxyphenyl) ethanol (natural antioxidant)

2 (4 ヒドロキシフェニル) エタノール (天然抗酸化剤)

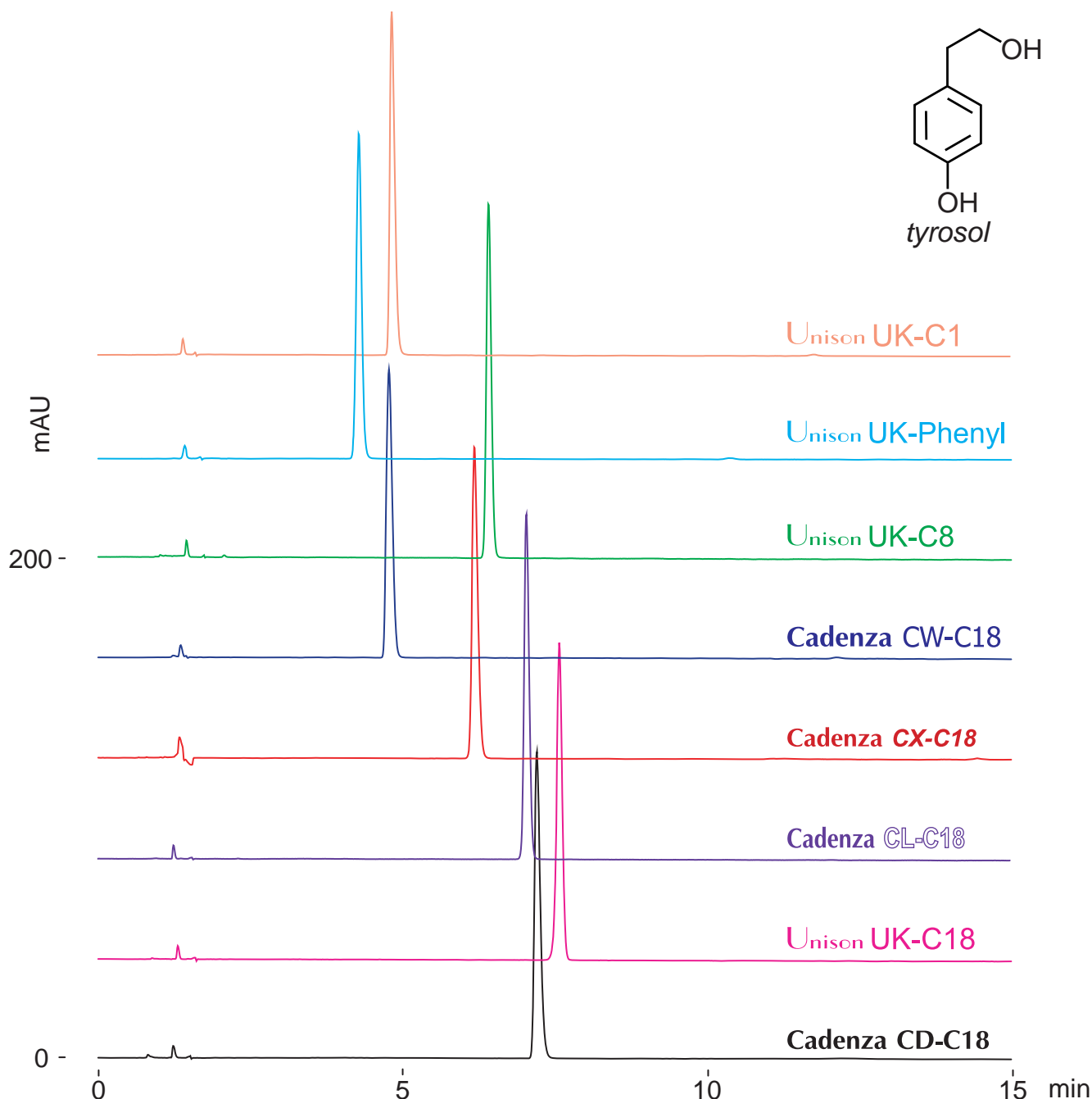

100 x 4.6 mm

A: water / acetonitrile / formic acid = 95 / 5 / 0.1

B: water / acetonitrile / formic acid = 50 / 50 / 0.1

0%B (0-1min), 0-30%B (1-15min)

1.0mL/min (9-11MPa), 37°C, 260nm

10µL (0.3mg/mL, water / formic acid = 100 / 1)