



Unison UK-C1
Unison UK-Phenyl
Unison UK-C8
Unison UK-C18
Cadenza CX-C18
Cadenza CW-C18
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

2- (4-Hydroxyphenyl) ethanol (natural antioxidant)

2 - (4 - ヒドロキシフェニル) エタノール (天然抗酸化剤)

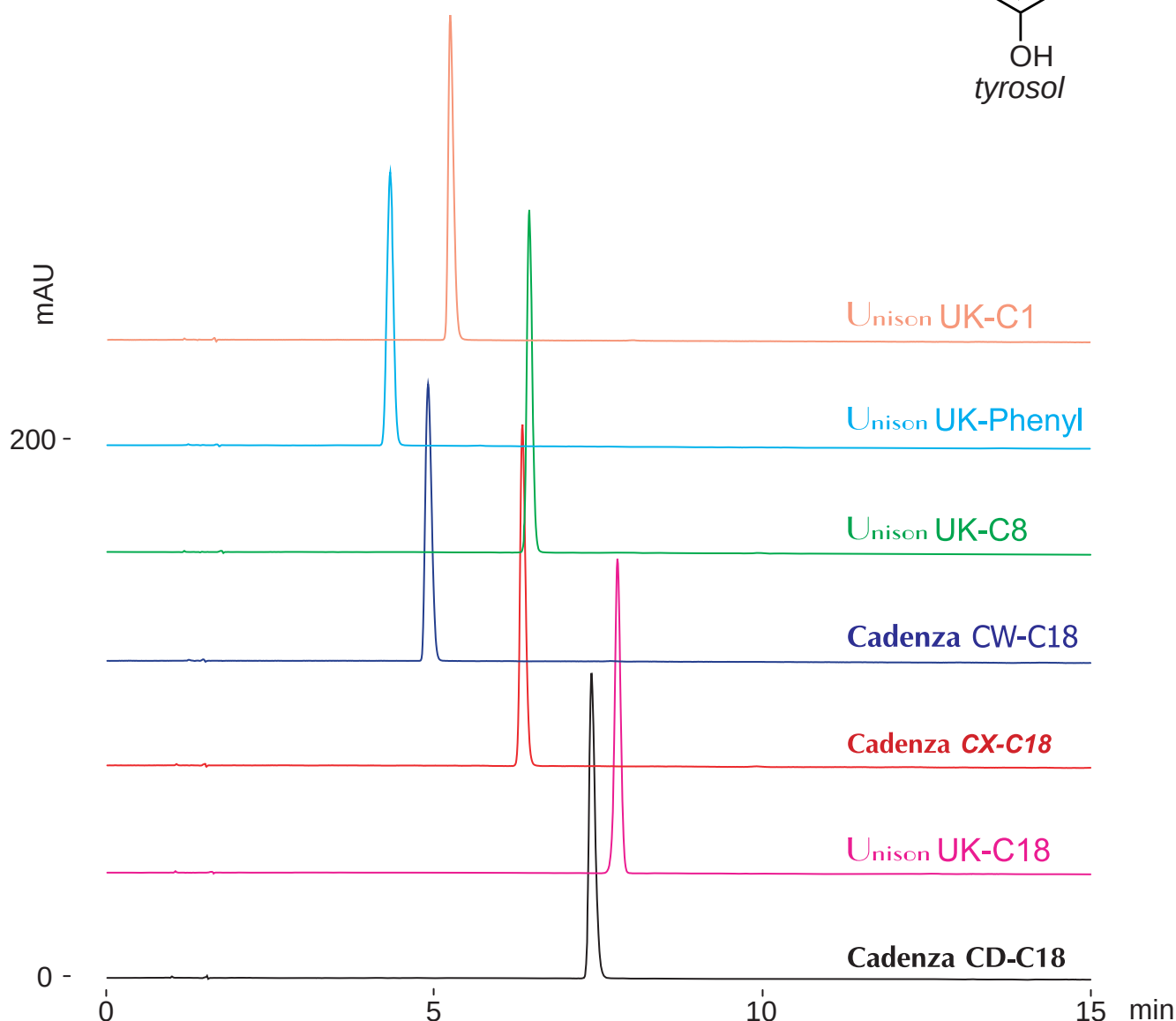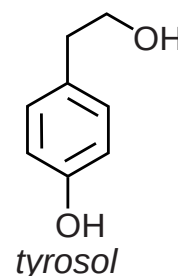

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 95 / 5

B: 100mM ammonium formate / acetonitrile = 50 / 50

0%B (0-1min), 0-30%B (1-15min)

1.0mL/min (9-11MPa), 37°C, 260nm

10μL (0.5mg/mL, 100mM ammonium formate)