
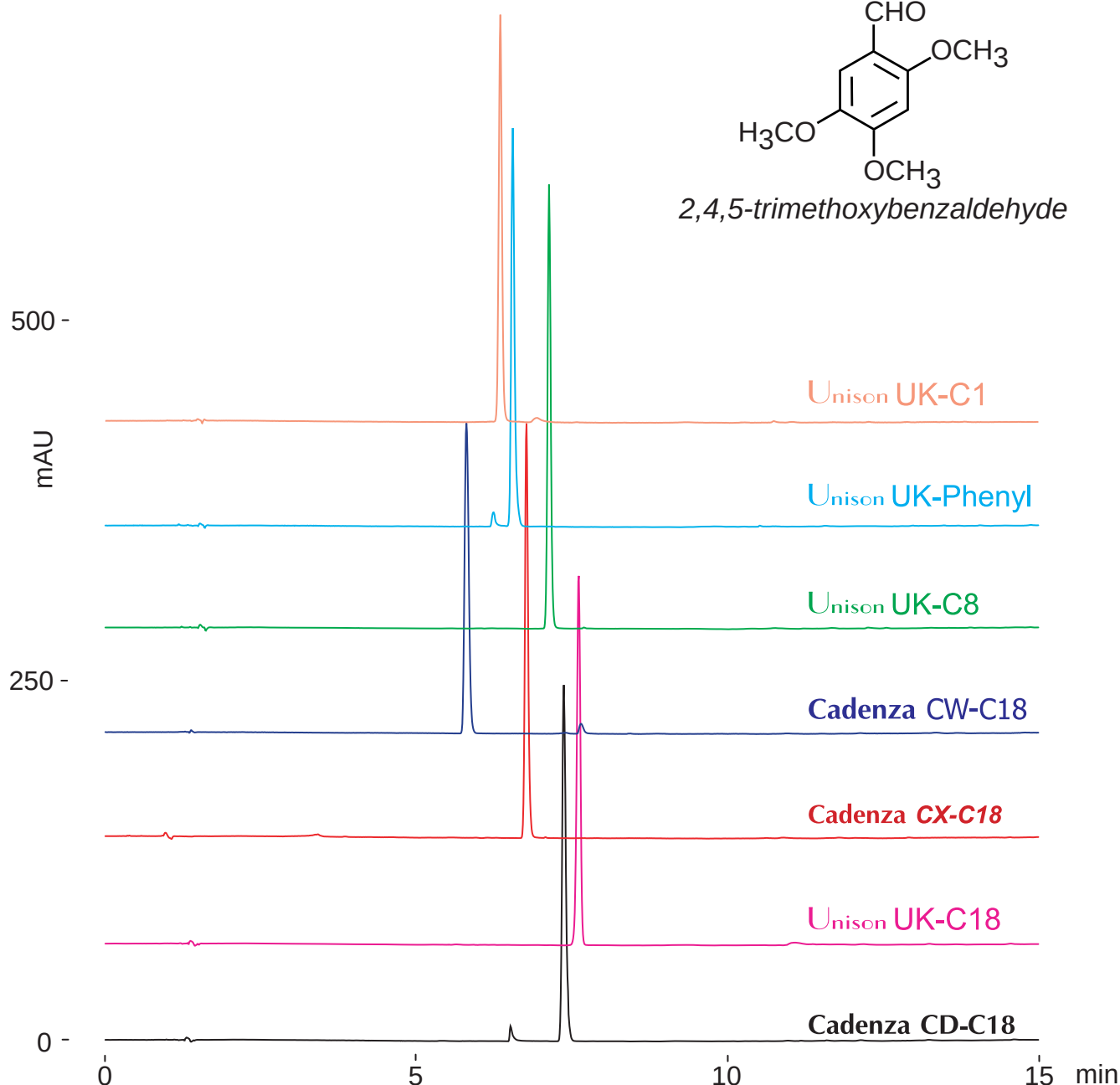


Unison UK-C1  
Unison UK-Phenyl  
Unison UK-C8  
Unison UK-C18  
Cadenza CX-C18  
Cadenza CW-C18  
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## 2,4,5-Trimethoxybenzaldehyde (Natural antibacterials) 2, 4, 5 -トリメトキシベンズアルデヒド (天然抗菌成分)

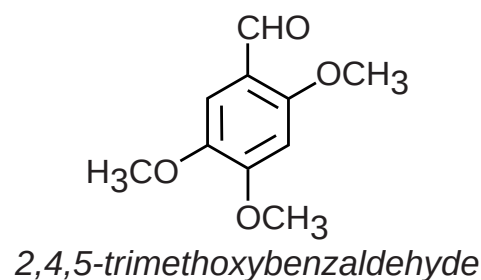

100 x 4.6 mm

A: water / acetonitrile / acetic acid = 90 / 10 / 0.5

B: water / acetonitrile / acetic acid = 10 / 90 / 0.5

0-80%B (0-15min)

1.0mL/min (9-11MPa), 37°C, 260nm

10μL (0.1mg/mL, water / acetonitrile / acetic acid = 80 / 20 / 1)