


**Cadenza CX-C18**

150 x 3 mm


**HPLC COLUMN APPLICATION**

(R) -Bis [3,5-bis (trifluoromethyl) phenyl] (pyrrolidin-2-yl) methanol  
(Proline type - organocatalyst)

(R) - ビス [3, 5 ビス (トリフルオロメチル) フェニル]  
(ピロリジン -2-イル) メタノール  
(プロリン型有機触媒)

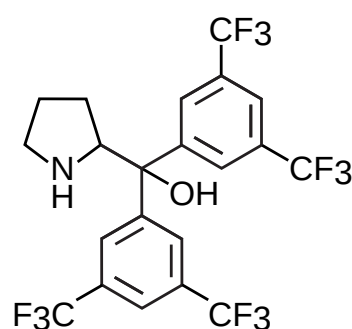

(R) -bis[ 3,5-bis (trifluoromethyl)  
phenyl] (pyrrolidin-2-yl) methanol

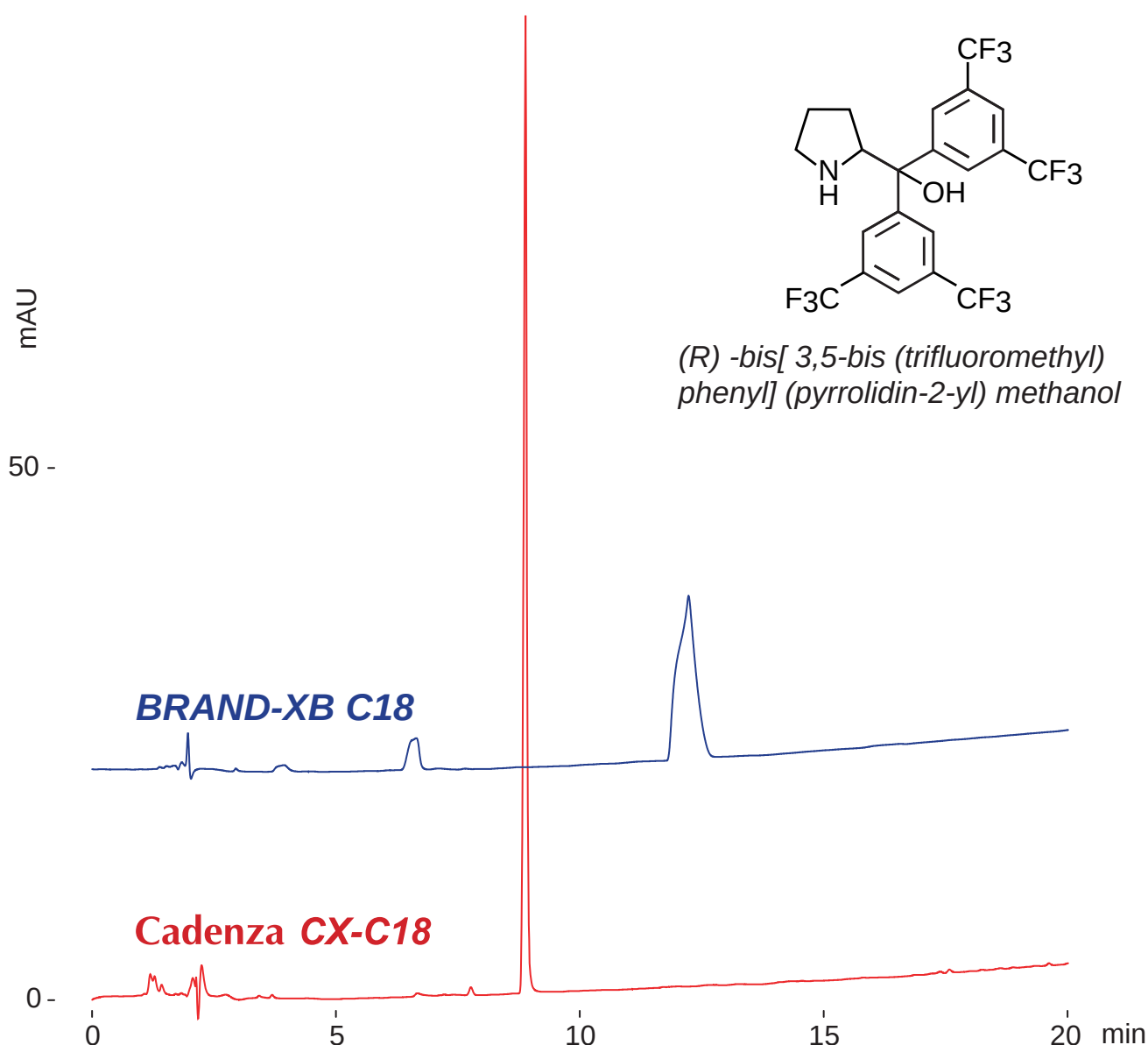

150 x 3 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.5

B: water / acetonitrile / formic acid = 10 / 90 / 0.5

10-100%B (0-15min)

0.4mL/min (11MPa), 37°C, 260nm

5uL (1.0mg/mL, water / acetonitrile / formic acid = 90 / 10 / 1)