



Unison UK-C1
Unison UK-Phenyl
Unison UK-C8
Unison UK-C18
Cadenza CX-C18
Cadenza CW-C18
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Isoliquiritigenin (Chalcones)

イソリキリチゲニン (カルコン類)

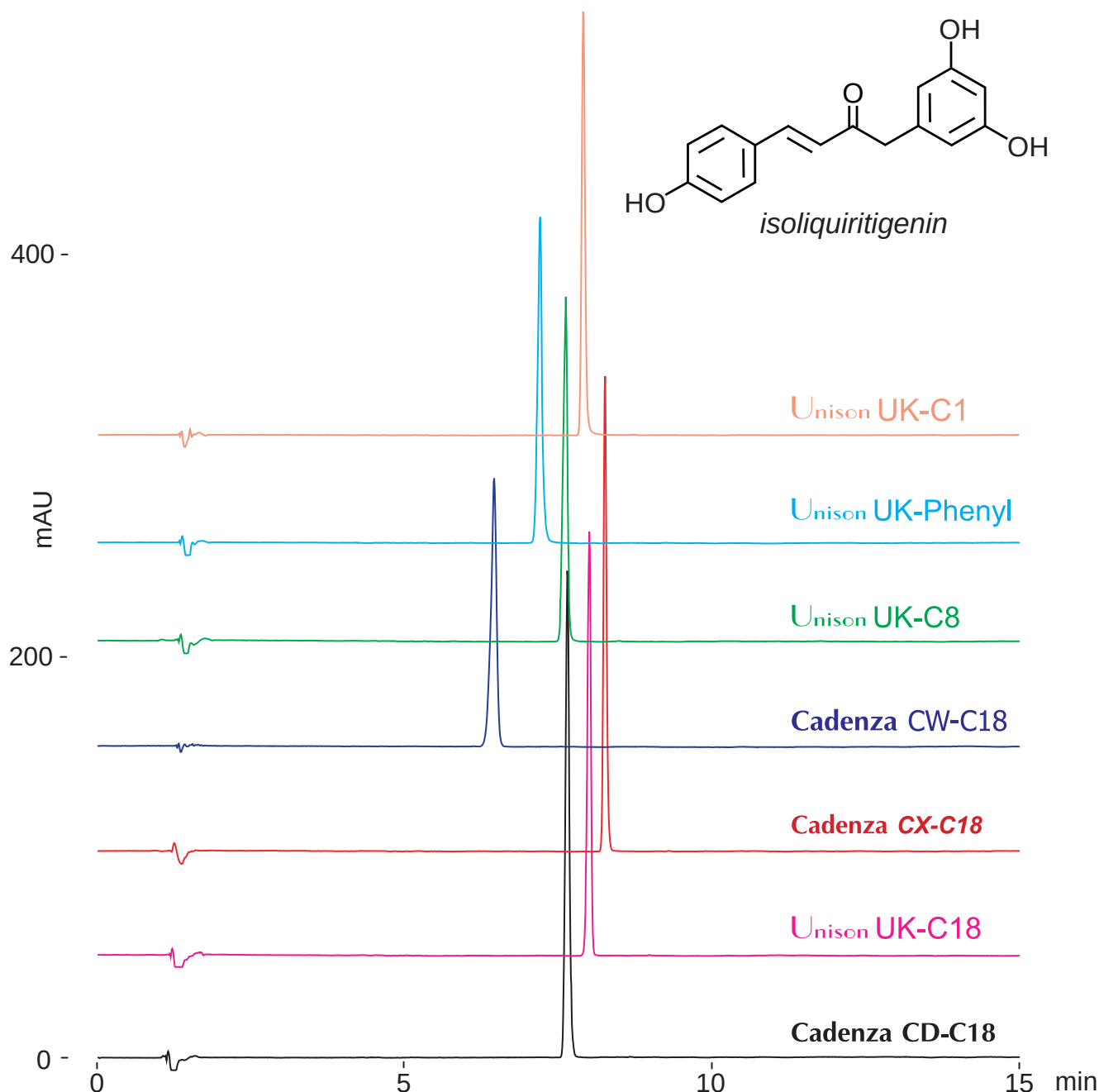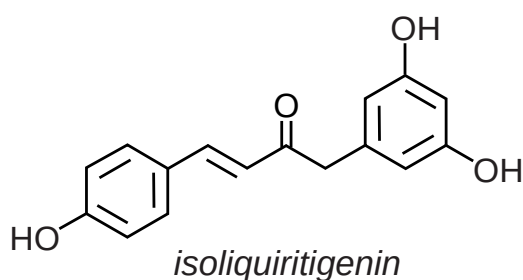

100 x 4.6 mm

A: water / acetonitrile / acetic acid = 90 / 10 / 0.5

B: water / acetonitrile / acetic acid = 10 / 90 / 0.5

5-100%B (0-15min)

1.0mL/min (9-12MPa), 37°C, 370nm

10μL (0.025mg/mL, water / ethanol / acetic acid = 10 / 90 / 1)