


Unison UK-Amino
Nardis ND-NX

150 x 4.6 mm

HPLC COLUMN APPLICATION

Iminoctadine(pesticides)

イミノクタジン (殺菌剤)

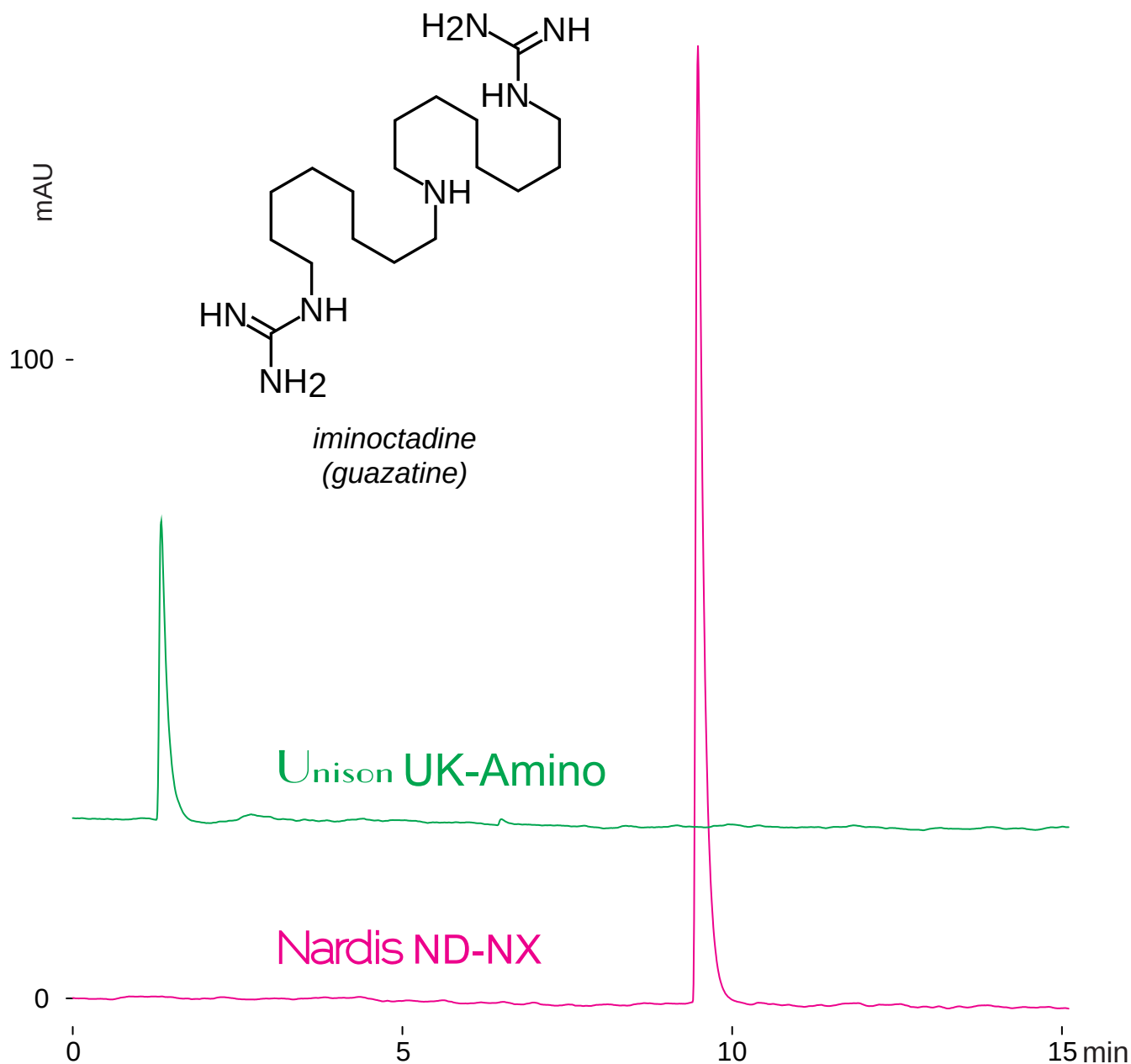

150 x 4.6 mm

A: 50mM ammonium formate / acetonitrile = 50 / 50

B: water / formic acid = 100 / 0.5

0-100%B (0-15min)

1.0mL/min (11-13MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10μL (1mg/mL, water)