
Cadenza *CX-C18* 150 x 3 mm **HPLC COLUMN APPLICATION**

(R)-(-)-2-[Hydroxy (diphenyl) methyl]-1-methylpyrrolidine (Proline-organocatalyst)

(R)-(-)-2-[ヒドロキシ (ジフェニル)メチル]-1 -メチルピロリジン (プロリン型有機触媒)

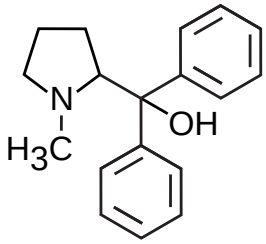

(R)-(-)-2-[hydroxy (diphenyl) methyl] -1-methylpyrrolidine

150 x 3 mm
A: water / formic acid = 100 / 0.5
B: water / acetonitrile / formic acid = 50 / 50 / 0.5
0-75%B (0-20min)
0.4mL/min (11MPa), 37°C, 260nm
2µL (2.0mg/mL, water / formic acid = 100 / 1)