
Dacapo DX-C18 100 x 4.6 mm HPLC COLUMN APPLICATION

(2S)-N-[(1S)-1-(Hydroxydiphenylmethyl)-3-methylbutyl] -2-pyrrolidinecarboxamide (Proline-organocatalyst)

(2S)-N-[(1S)-1-(ヒドロキシジフェニルメチル)-3-メチルブチル] -2-ピロリジンカルボキサミド (プロリン型有機触媒)

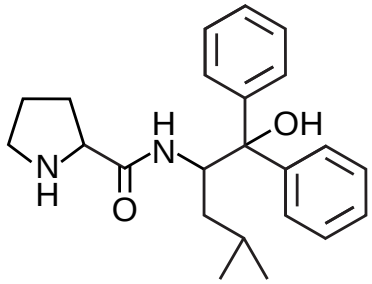

(2S)-N-[(1S)-1-(hydroxydiphenylmethyl) -3-methylbutyl]-2-pyrrolidinecarboxamide

Dacapo DX-C18, 100 x 4.6 mm
A: 10mM NH4OH / acetonitrile = 90/ 10
B: 10mM NH4OH / acetonitrile = 10 / 90
5-100%B (0-15min)
1mL/min (18MPa),37°C, 260nm
10µL (2.0mg/mL, 100mM ammonium formate / acetonitrile = 67 / 33)