
Cadenza CD-C18 100 x 3 mm **HPLC COLUMN APPLICATION**

LC-MS: beta-Carboline alkaloids

β -カルボリンアルカロイド類の LC-MS分析

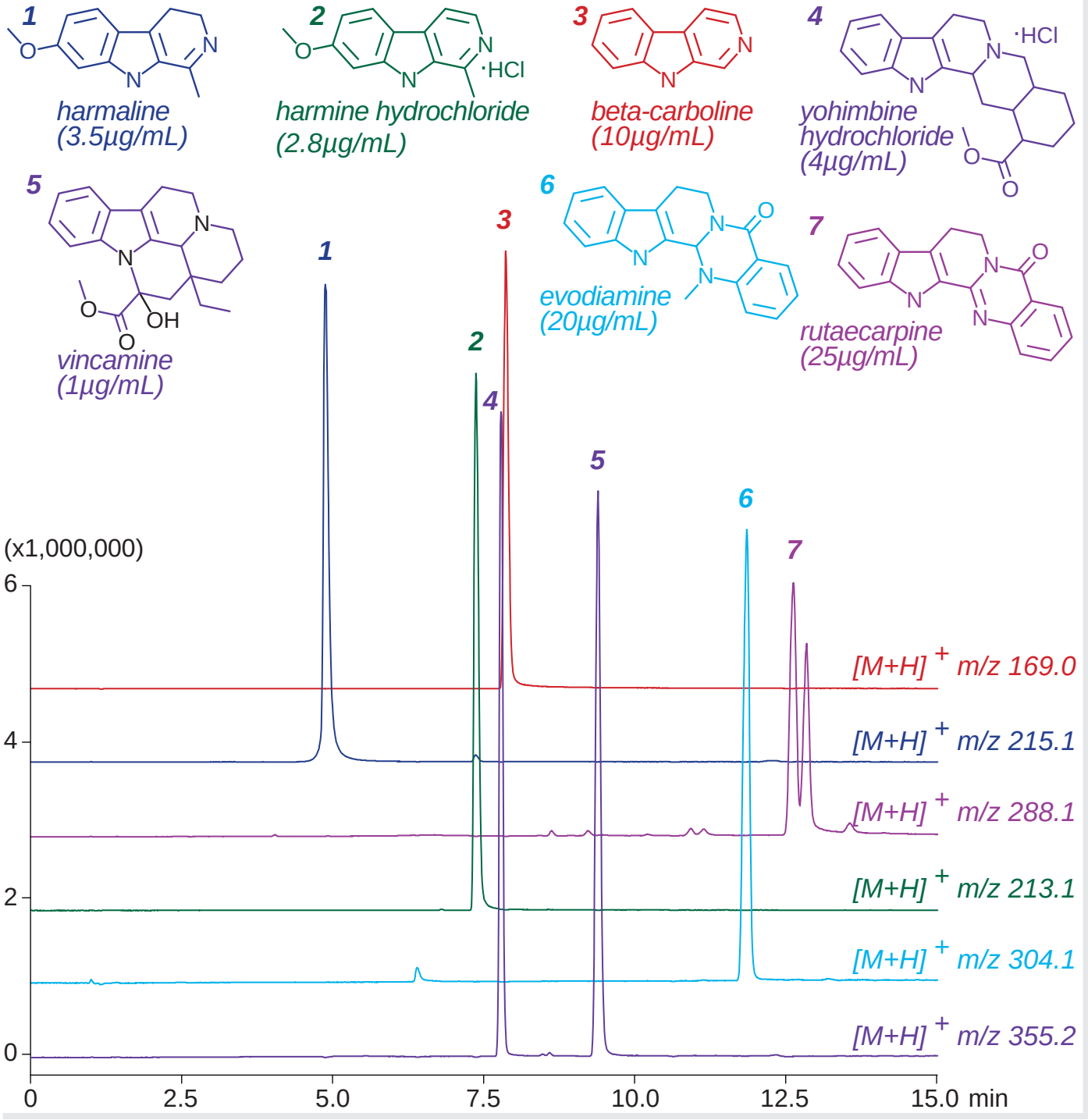

Cadenza CD-C18, 100 x 4.6 mm
A: acetonitrile / 100mM ammonium acetate = 10 / 90
B: acetonitrile / 100mM ammonium acetate = 95 / 5
5 -80%B (0-15min), 100%B (15-17min)
0.4mL/min (9MPa), 37°C, Single Quad. MS (ESI)
2µL (1-25µg/mL, ACN / 100mM AcONH4 = 50 / 50)