
Scherzo SS-C18
Scherzo SM-C18
Scherzo SW-C18

100 x 4.6 mm HPLC COLUMN APPLICATION

(S)-5-(Pyrrolidin-2-yl)-1H-tetrazole (Proline-type organocatalyst)

(S)-5-(ピロリジン-2-イル)-1H-テトラゾール (プロリン型有機触媒)

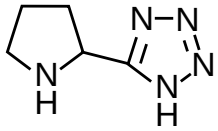

(S)-5-(pyrrolidin-2-yl)-1H-tetrazole

100 x 4.6 mm
A: water / acetonitrile / formic acid = 95 / 5 / 0.1
B: 100mM ammonium formate / acetonitrile = 50 / 50
1.0mL/min (9MPa), 37°C, ELSD (spray chamber 30°C, drift tube 60°C)
10µL (2.0mg/mL, water / formic acid = 100 / 1)