

Scherzo SM-C18
Nardis ND-RX 100 x 3 mm HPLC COLUMN APPLICATION

LC-MS: Chloroquine analogues (Antimalarial)


クロロキン類縁体（マラリア薬）の LC-MS分析

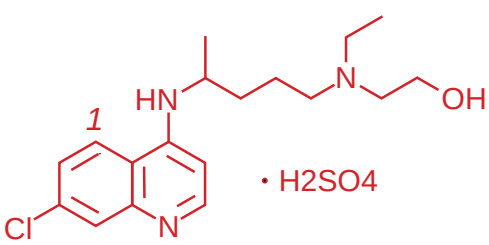

hydroxychloroquine sulfate

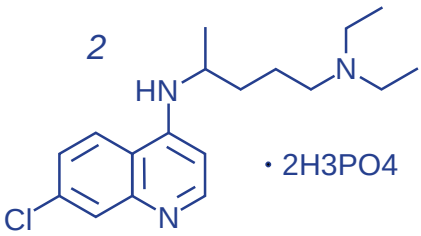

chloroquine diphosphate

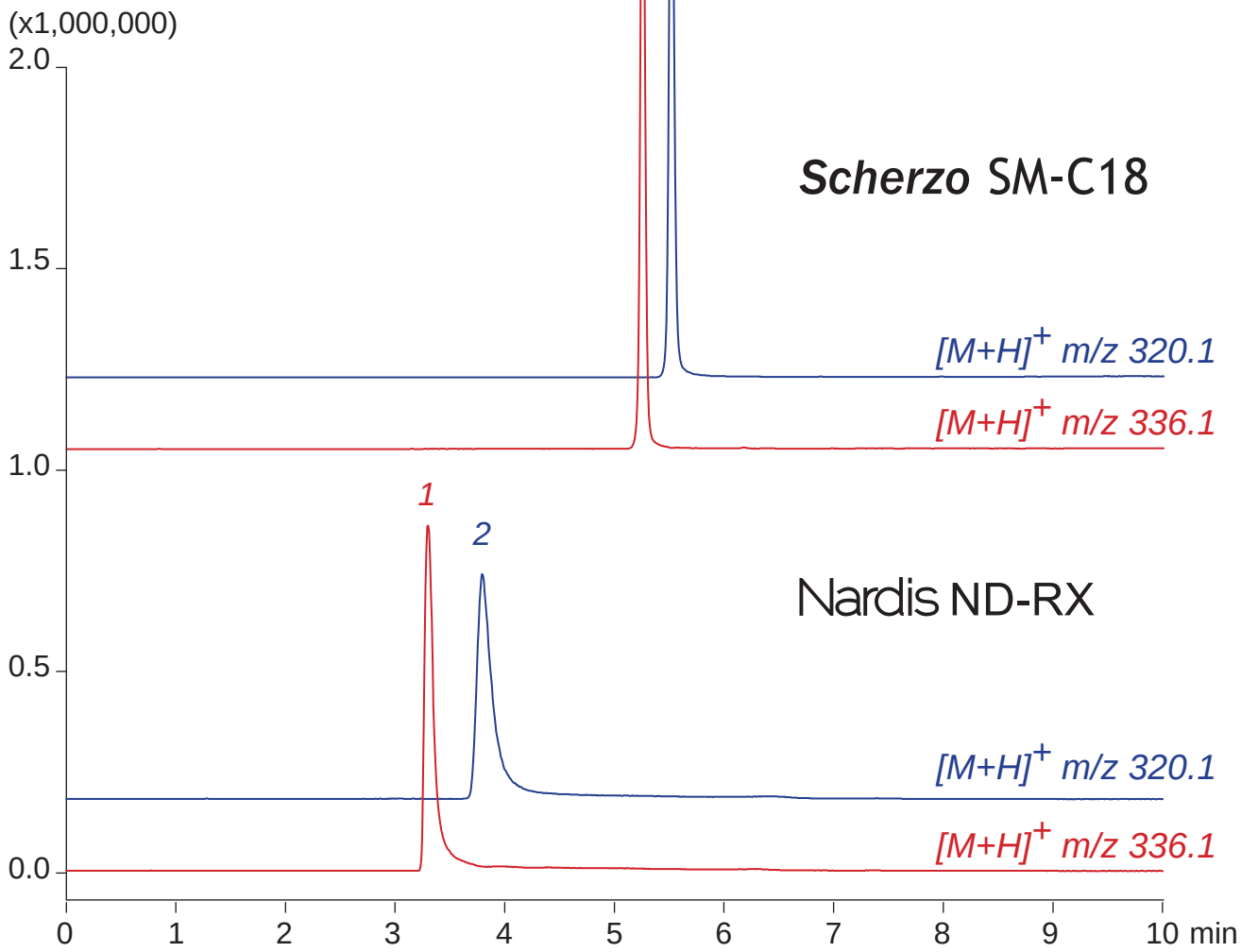

100 x 3 mm
A: water / formic acid = 100 / 0.1
B: 100mM ammonium formate / acetonitrile = 50 / 50
0-100%B (0-15min)
0.4mL/min (6-11MPa), 37°C, Single Quad. MS (ESI)
5µL (1µg/mL, water / formic acid = 100 / 1)