

Unison UK-Phenyl

Nardis ND-RX


Nardis NF-PFP/Biphenyl

150 x 4.6 mm

HPLC COLUMN APPLICATION

## Chloromandelic acid isomers

クロロマンデル酸異性体

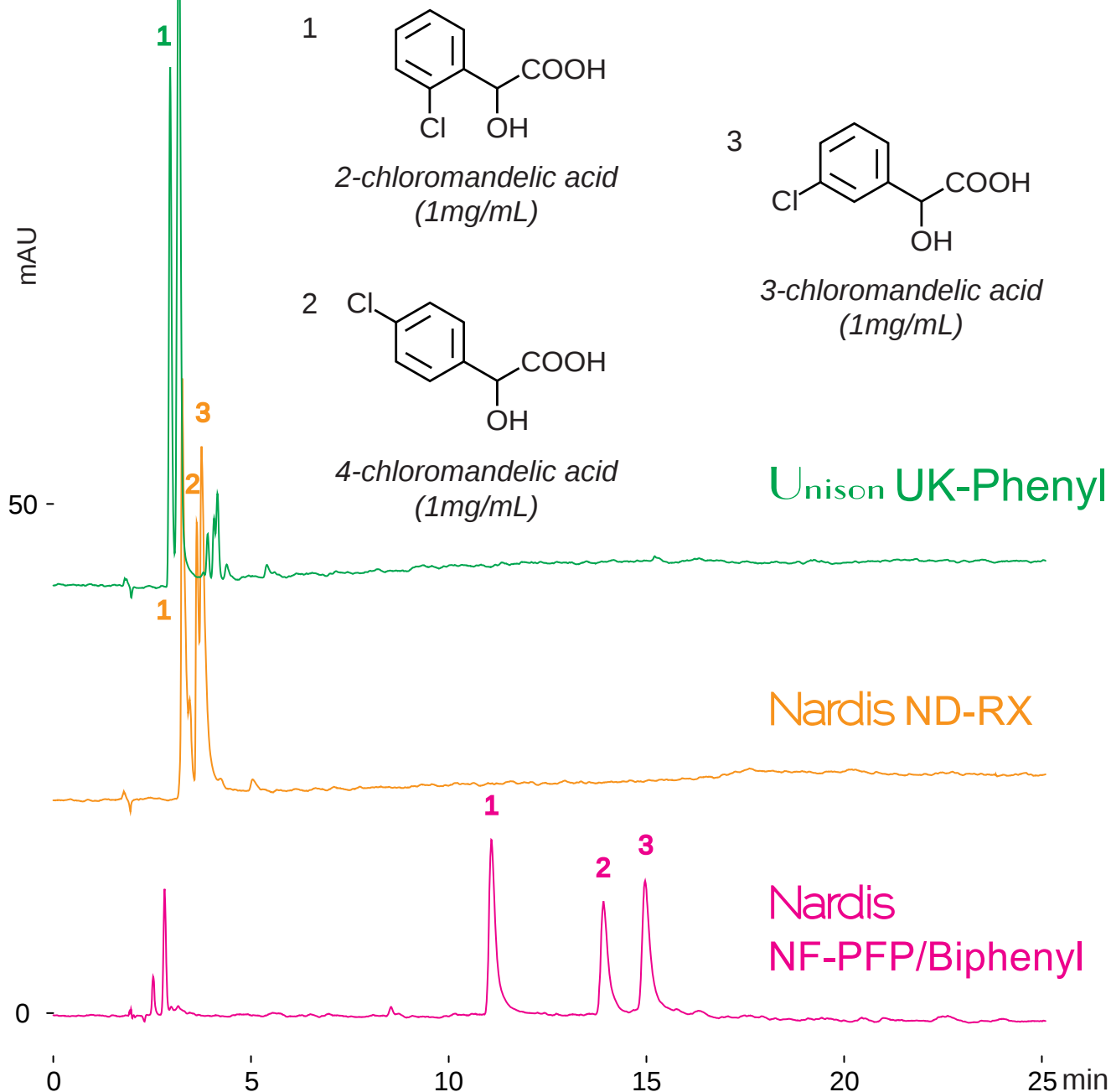

150 x 4.6 mm

A: water / methanol / formic acid = 50 / 50 / 0.5

B: methanol / formic acid = 100 / 0.5

0-100%B (0-15min), 100%B (15-25min)

1.0mL/min (15-17MPa), 37°C, 270nm

10μL (1mg/mL, methanol)