

Nardis NF-PFP/Biphenyl

Nardis ND-RX



Unison UK-Phenyl

150 x 4.6 mm


HPLC COLUMN APPLICATION

## Nonsteroidal Antiinflammatory Drugs (NSAIDs)


非ステロイド性抗炎症薬

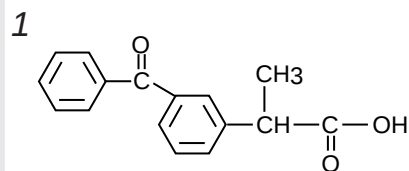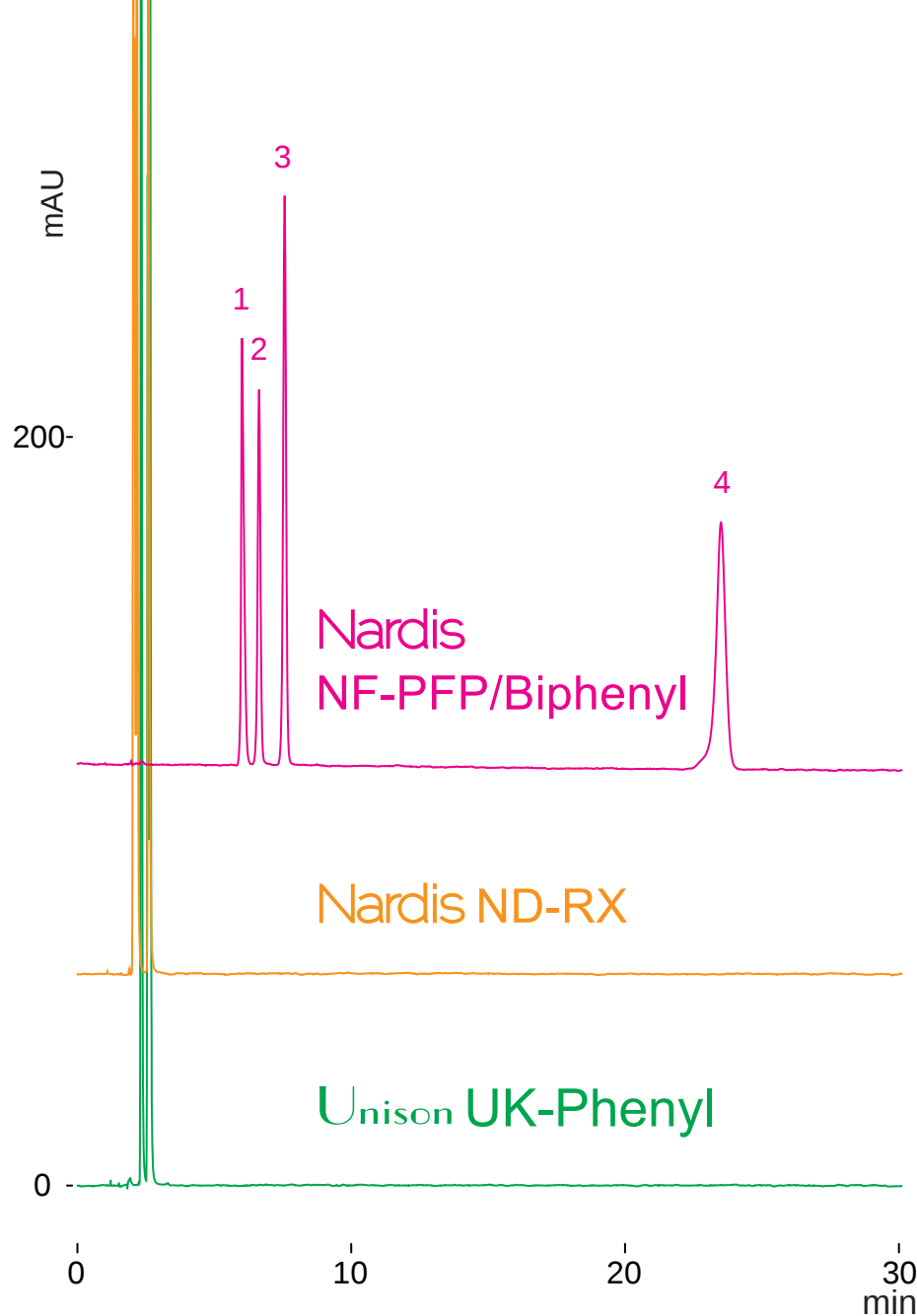


*ketoprofen*  
(0.08mg/mL)

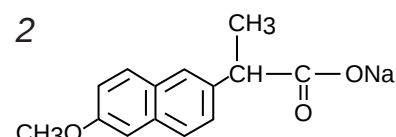

*naproxen sodium*  
(0.15mg/mL)

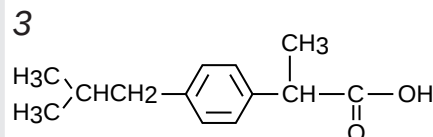

*ibuprofen*  
(3mg/mL)

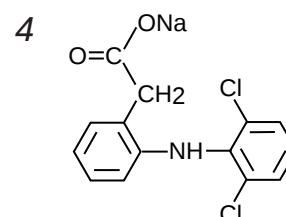

*diclofenac sodium*  
(0.3mg/mL)

150 x 4.6 mm

A: water / acetonitrile / acetic acid = 50 / 50 / 0.5

B: water / acetonitrile / acetic acid = 10 / 90 / 0.5

50-100%B (0-15min), 100%B (15-30min)

1.0mL/min (7-10MPa), 37°C, 260nm

10μL (3-0.08mg/mL, ACN/MeOH/water/AcOH=70/20/10/2)