


Presto FF-C18

100 x 3 mm

HPLC COLUMN APPLICATION

Non-Porous ODS-MS:human Insulin (aggregationsolid matter)

変性ヒトインスリン固形物のノンポーラスODSカラム-MS分析

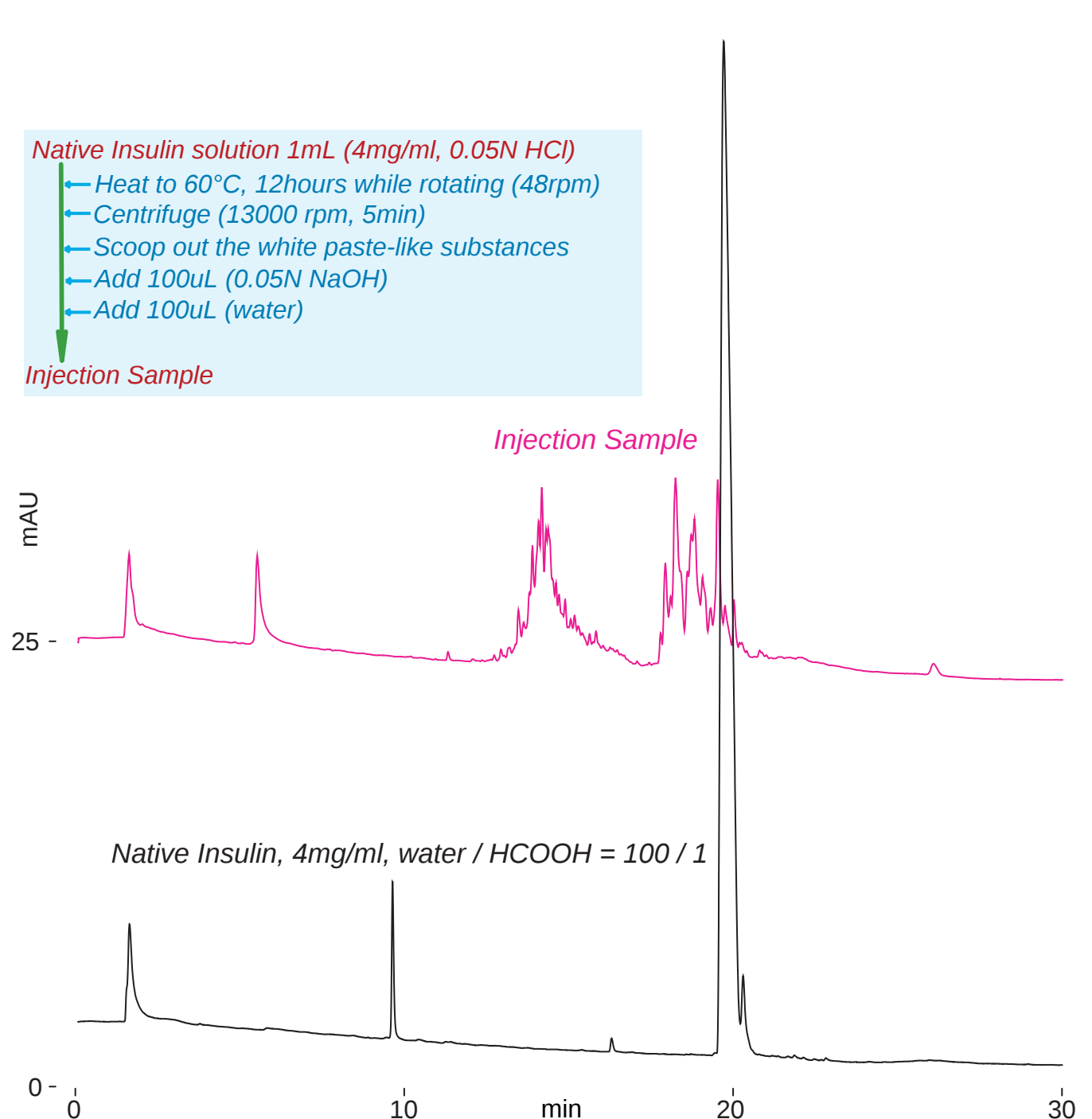

Presto FF-C18, 100 x 3mm

A: water / trifluoroacetic acid = 100 / 0.1

B: acetonitrile / trifluoroacetic acid = 100 / 0.1

5-50%B (0-30min)

0.2mL/min (15.7MPa), 37°C, 280nm

5uL