


Presto FF-C18

75 x 3 mm

HPLC COLUMN APPLICATION

Non-Porous ODS-MS:rHSA

遺伝子組み換えヒト血清アルブミン (rHSA) のノンポーラスODSカラム-MS分析

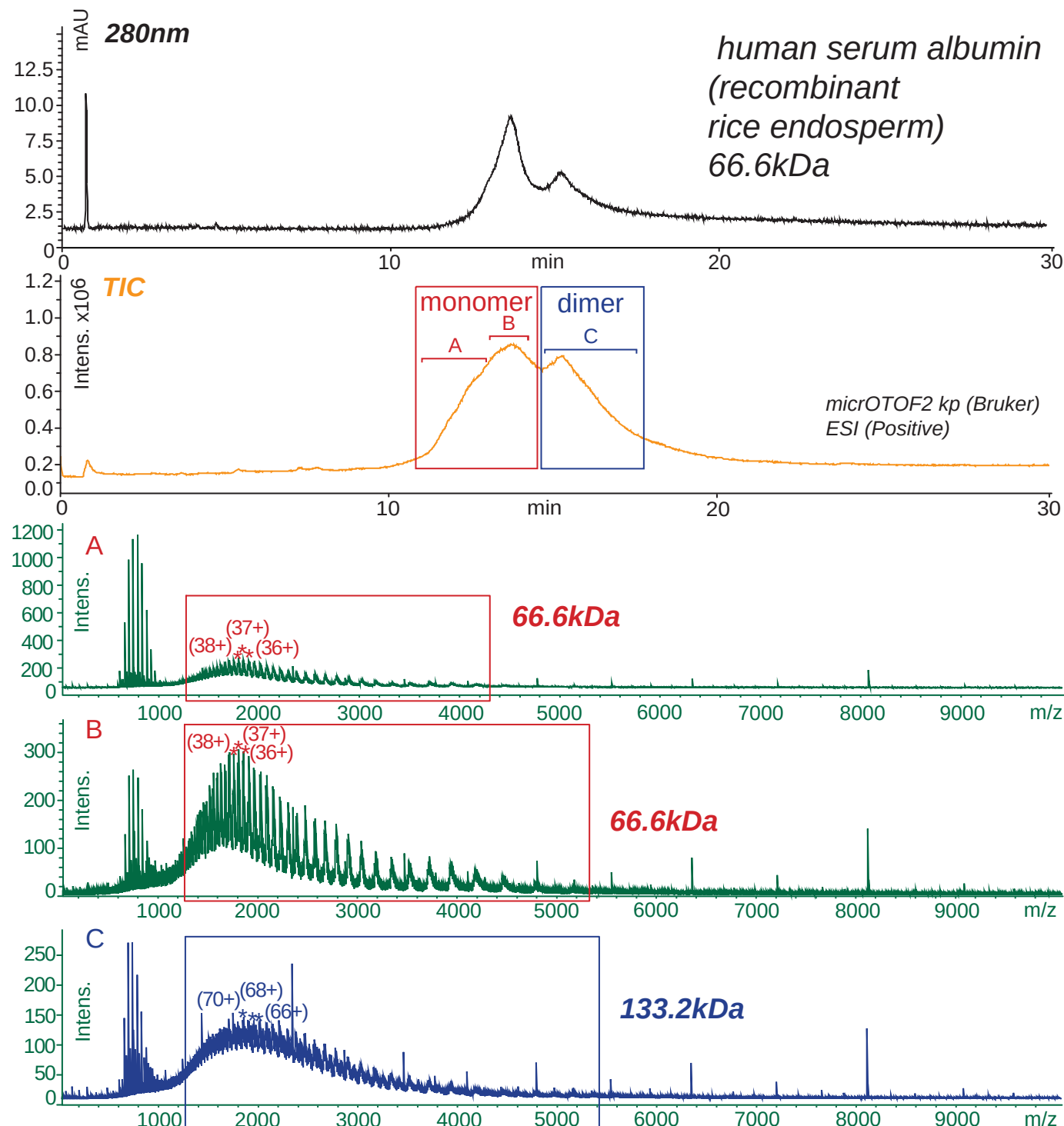

Presto FF-C18, 75 x 3mm

A: 100mM ammonium formate / acetonitrile = 95 / 5

B: 100mM ammonium formate / acetonitrile = 20 / 80

27-62%B (0-30min)

0.3mL/min (20MPa), 37°C

5uL (5mg/mL, 100mM ammonium formate)