


Presto FF-C18

75 x 3 mm

HPLC COLUMN APPLICATION

Non-Porous ODS-MS:Papain パパインのノンポーラスODSカラム - MS分析

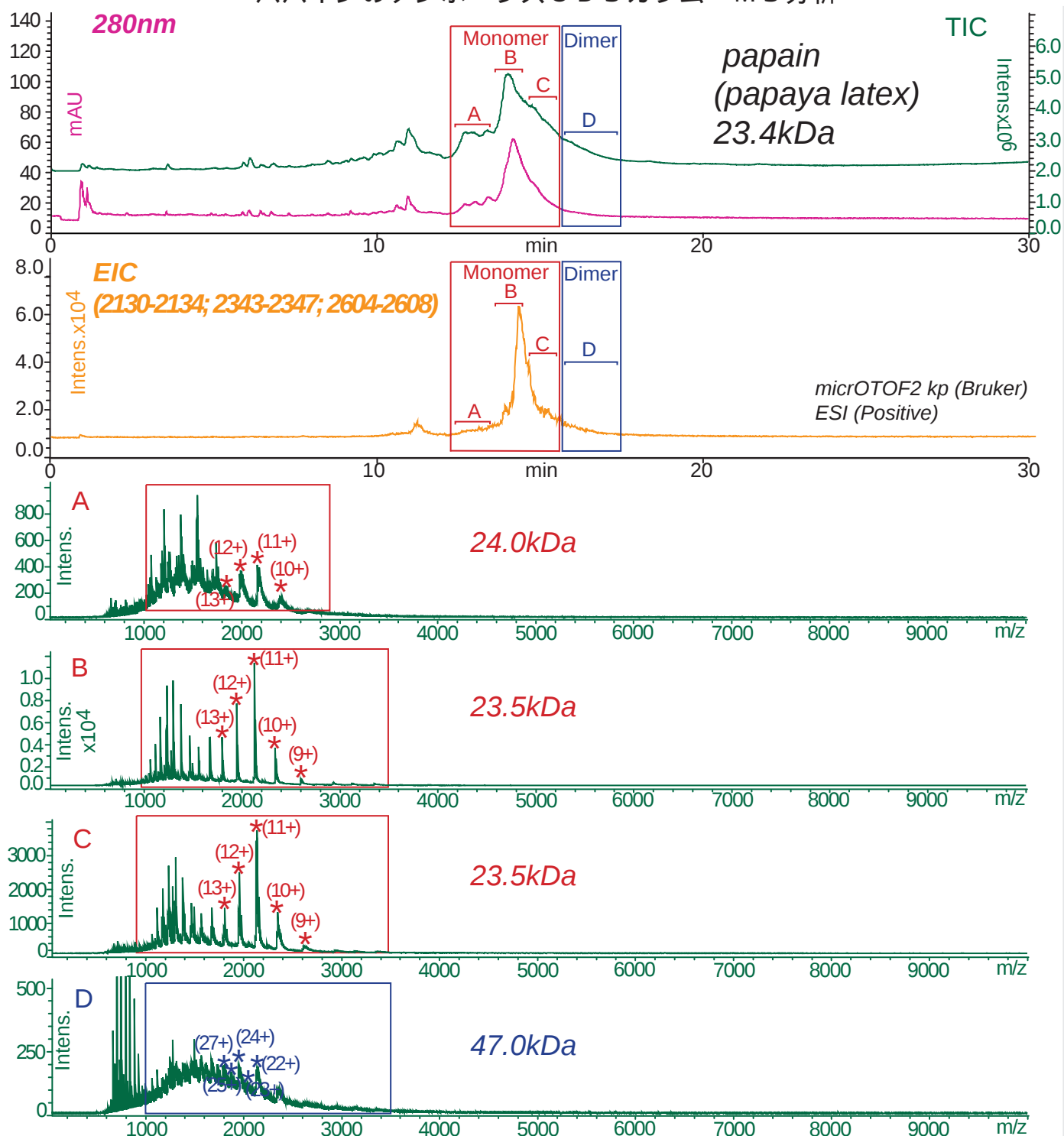

Presto FF-C18, 75 x 3mm

A: water / acetonitrile / formic acid = 95 / 5 / 0.5

B: water / acetonitrile / formic acid = 20 / 80 / 0.5

0-80%B (0-30min)

0.2mL/min (21MPa), 37°C

5 μ L (5mg/mL, water / formic acid = 100 / 1)