


Presto FF-C18

75 x 3 mm

HPLC COLUMN APPLICATION

Non-Porous ODS-MS:alpha-chymotrypsinogen A from bovine pancreas

α - キモトリプシノーゲン A (ウシ膵臓由来) のノンポーラス ODS カラム - MS 分析

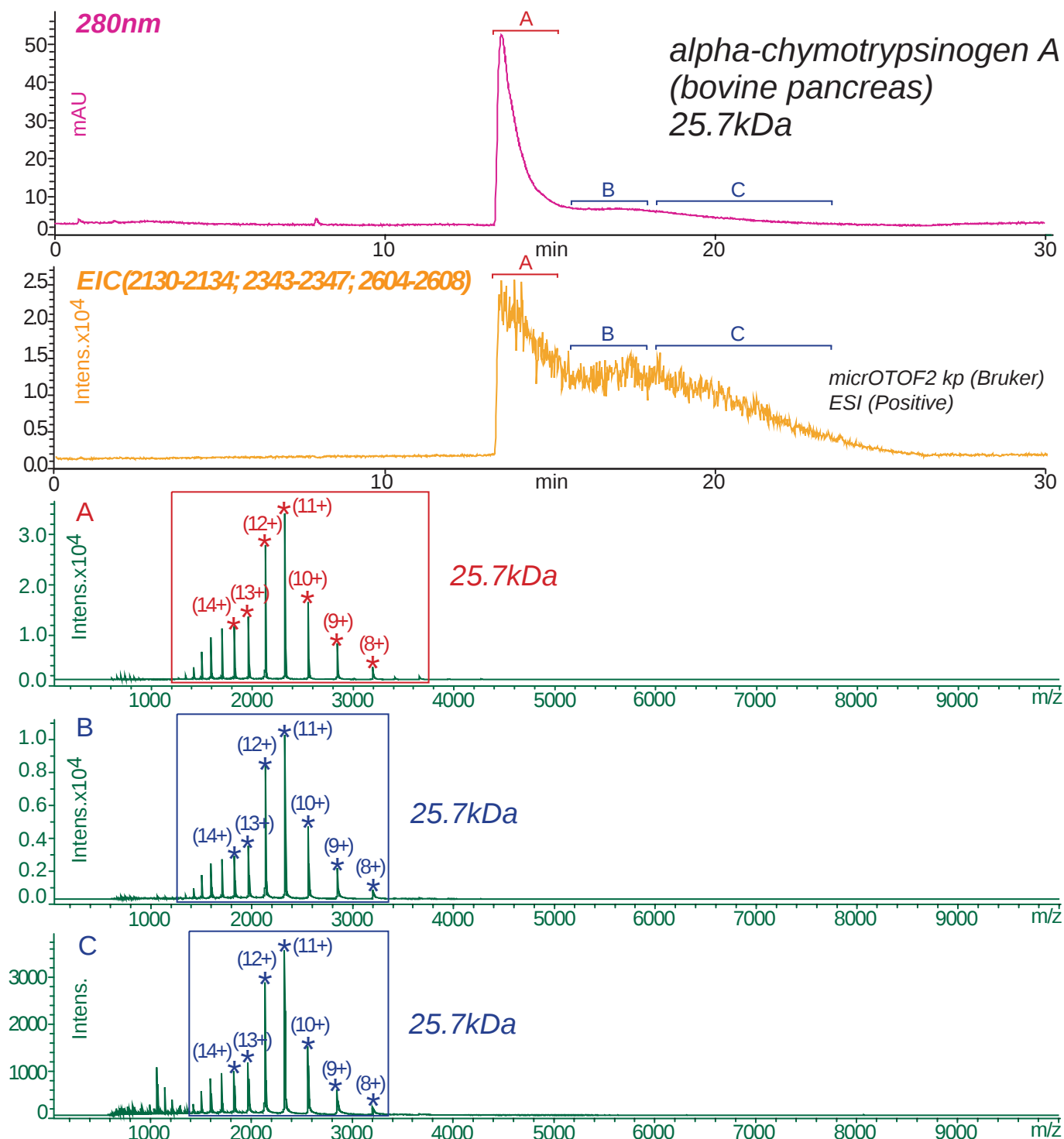

Presto FF-C18, 75 x 3mm

A: water / acetonitrile / formic acid = 95 / 5 / 0.5

B: water / acetonitrile / formic acid = 20 / 80 / 0.5

0-100%B (0-20min), 100%B (20-30min)

0.3mL/min (21MPa), 37°C

5μL (5mg/mL, water / formic acid = 100 / 1)