


Nardis NF-PFP/Biphenyl 250 x 4.6 mm


HPLC COLUMN APPLICATION

Methylacetophenone isomers


メチルアセトフェノン異性体

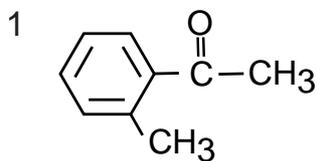


*o*-methylacetophenone  
(0.4mg/mL)

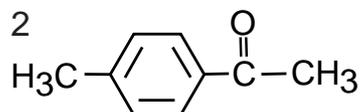

*p*-methylacetophenone  
(0.05mg/mL)

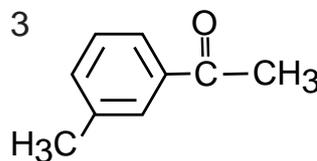

*m*-methylacetophenone  
(0.2mg/mL)

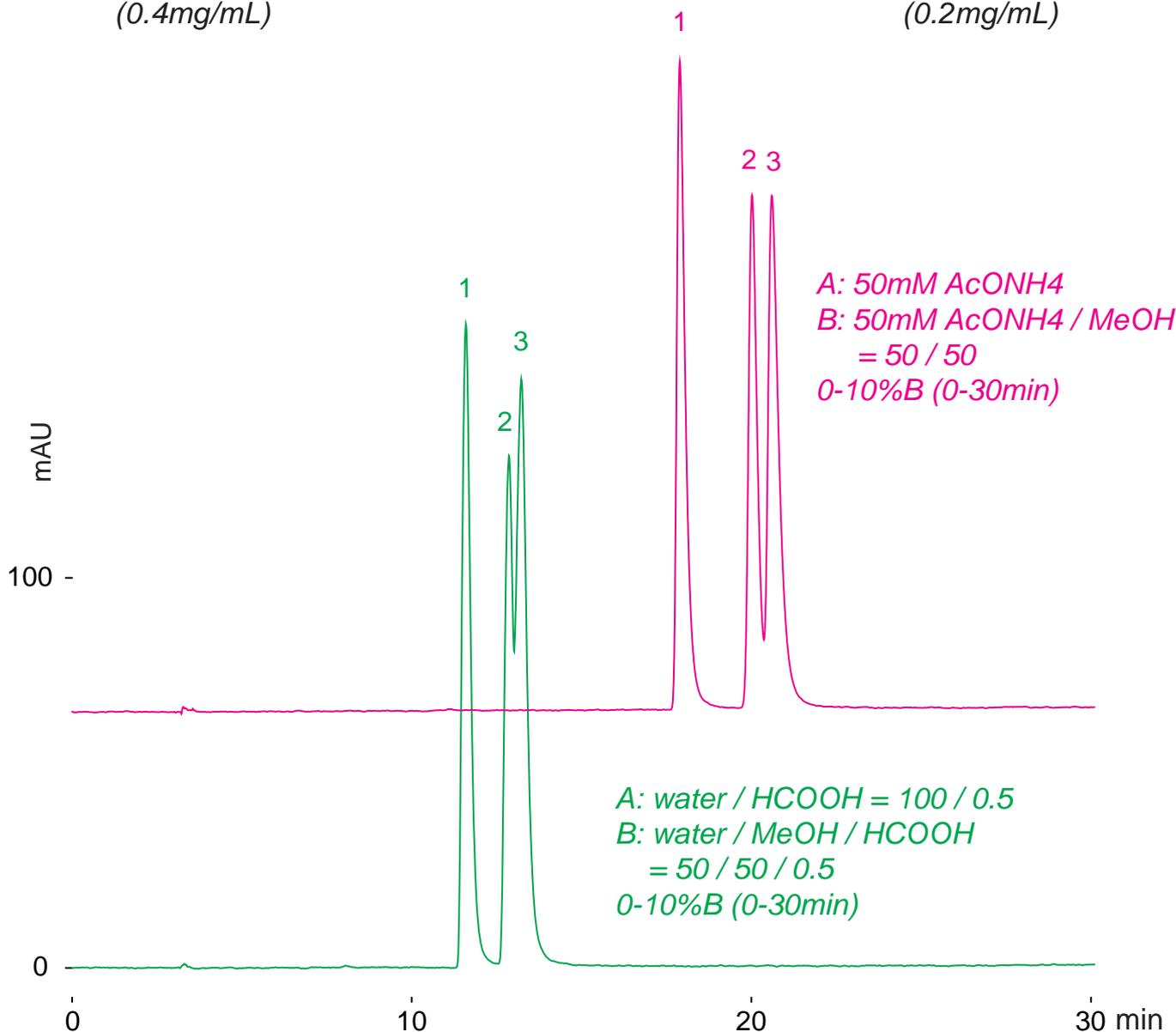

250 x 4.6 mm  
1.0mL/min (13-14MPa), 60°C, 260nm  
10µL (0.05-0.4mg/mL, methanol)