


Intrada Organic Acid

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: L-Threonic acid calcium salt (Supplement component)

L-トレオン酸 カルシウム塩 (サプリメント成分) の LC-MS 分析

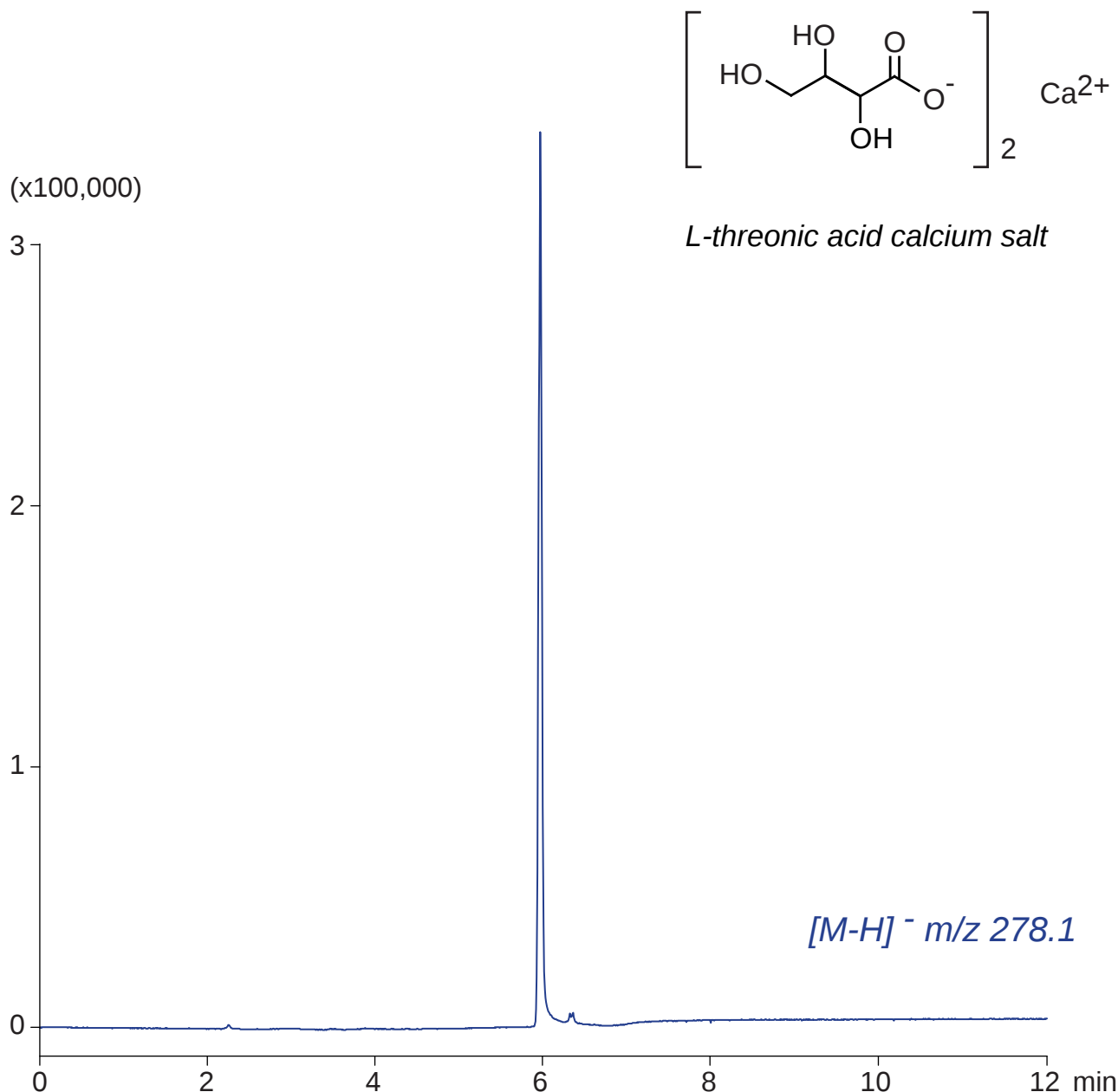

Intrada Organic Acid, 150 x 2 mm

A: acetonitrile / formic acid = 100 / 0.1

B: acetonitrile / 100mM ammonium formate = 10 / 90

10-100%B (0-6min), 100%B (6-12min)

0.2mL/min (5MPa), 37°C, Single Quad. MS (ESI)

5μL (100μg/mL, ACN / water / formic acid = 70 / 30 / 0.1)