




## 充てん剤バッチ間再現性

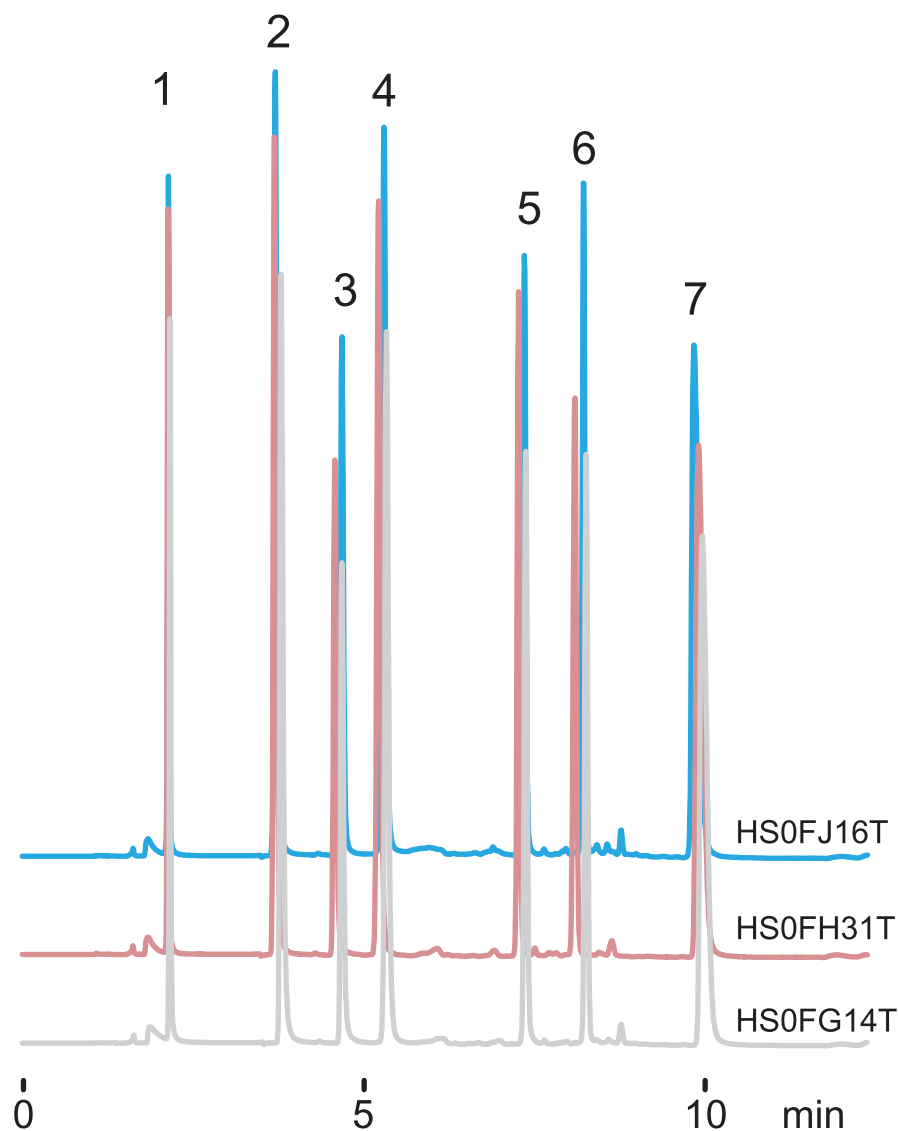

### Cadenza **HS-C18**

- 1 acyclovir (0.10ug)
- 2 theophylline (0.32ug)
- 3 ketoprofen (0.11ug)
- 4 slufamethazine (0.20ug)
- 5 carbamazepine (0.32ug)
- 6 progesterone (0.29ug)
- 7 nortriptyline (1.2ug)

Cadenza HS-C18, 150 x 3 mm

A: 10mM ammonium acetate

B: acetonitrile

0-75%B (0-10min)

0.4 mL/min

37 deg.C, 260 nm,

1.6 uL